



Alanine

Column: Chirex 3126 (D)-penicillamine**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-3126-E0**Elution Type:** Isocratic**Eluent A:** 1 mM Copper (II) sulfate in water**Pct. A:** 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

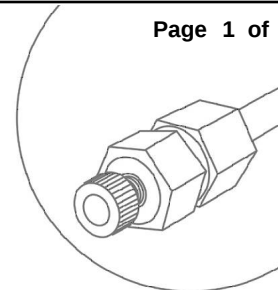

Products used in this Application:

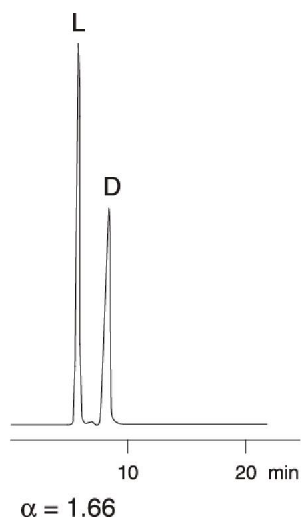CHIREX[®]**ANALYTES:**

- 1 L-Alanine
- 2 D-Alanine

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.