

PAH with Surrogate Standards on Envirosep-PP (UV), EPA Method 610

Page 1 of 1

Column: EnviroSep-PP

Dimensions: 125 x 4.6 mm ID

Order No: 00E-3029-E0

Elution Type: Gradient

Eluent A: Water



Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	25	0	100

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

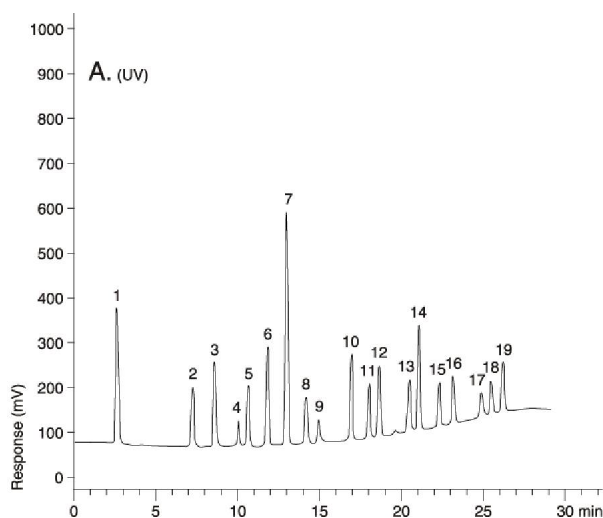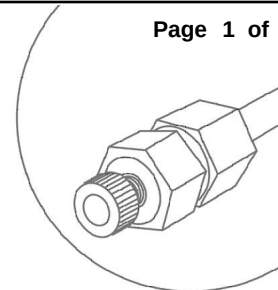

ANALYTES:

- | | |
|-------------------------|----------------------------|
| 1 Nitrobenzene | 15 Benzo(k)fluoranthene |
| 2 Naphthalene | 16 Benzo(a)pyrene |
| 3 Acenaphthylene | 17 Dibenzo(a,h)anthracene |
| 4 Acenaphthene | 18 Benzo(g,h,i)perylene |
| 5 Fluorene | 19 Indeno(1,2,3-c,d)pyrene |
| 6 Phenanthrene | |
| 7 Anthracene | |
| 8 Fluoranthene | |
| 9 Pyrene | |
| 10 p-Terphenyl | |
| 11 Benzo(a)anthracene | |
| 12 Chrysene | |
| 13 Benzo(e)pyrene | |
| 14 Benzo(b)fluoranthene | |

