

PAH with Surrogate Standards on Envirosep-PP (Flu), EPA Method 610

Page 1 of 1

Column: EnviroSep-PP

Dimensions: 125 x 4.6 mm ID

Order No: 00E-3029-E0

Elution Type: Gradient

Eluent A: Water



Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	25	0	100

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: Fluorescence @ 240 nm / 425 nm (ambient)

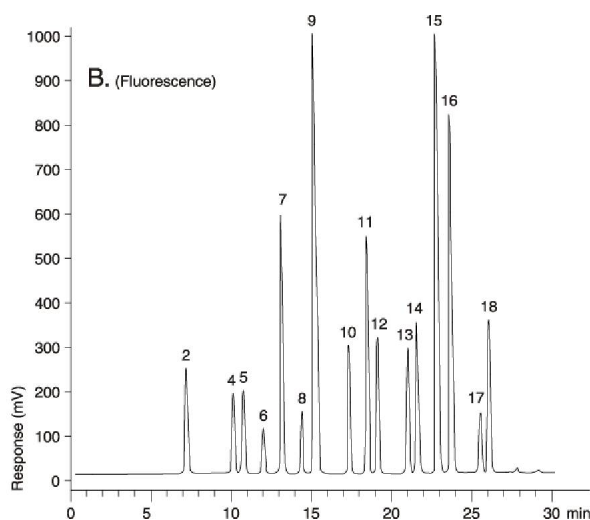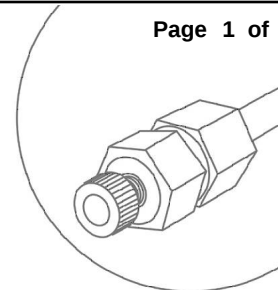

ANALYTES:

- | | |
|-------------------------|----------------------------|
| 1 Nitrobenzene | 15 Benzo(k)fluoranthene |
| 2 Naphthalene | 16 Benzo(a)pyrene |
| 3 Acenaphthylene | 17 Dibenzo(a,h)anthracene |
| 4 Acenaphthene | 18 Benzo(g,h,i)perylene |
| 5 Flourene | 19 Indeno(1,2,3-c,d)pyrene |
| 6 Phenanthrene | |
| 7 Anthracene | |
| 8 Fluoranthene | |
| 9 Pyrene | |
| 10 p-Terphenyl | |
| 11 Benzo(a)anthracene | |
| 12 Chrysene | |
| 13 Benzo(e)pyrene | |
| 14 Benzo(b)fluoranthene | |

©Copyright 2000-2006 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


 ...breaking with tradition[™]

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au