


Food Softeners

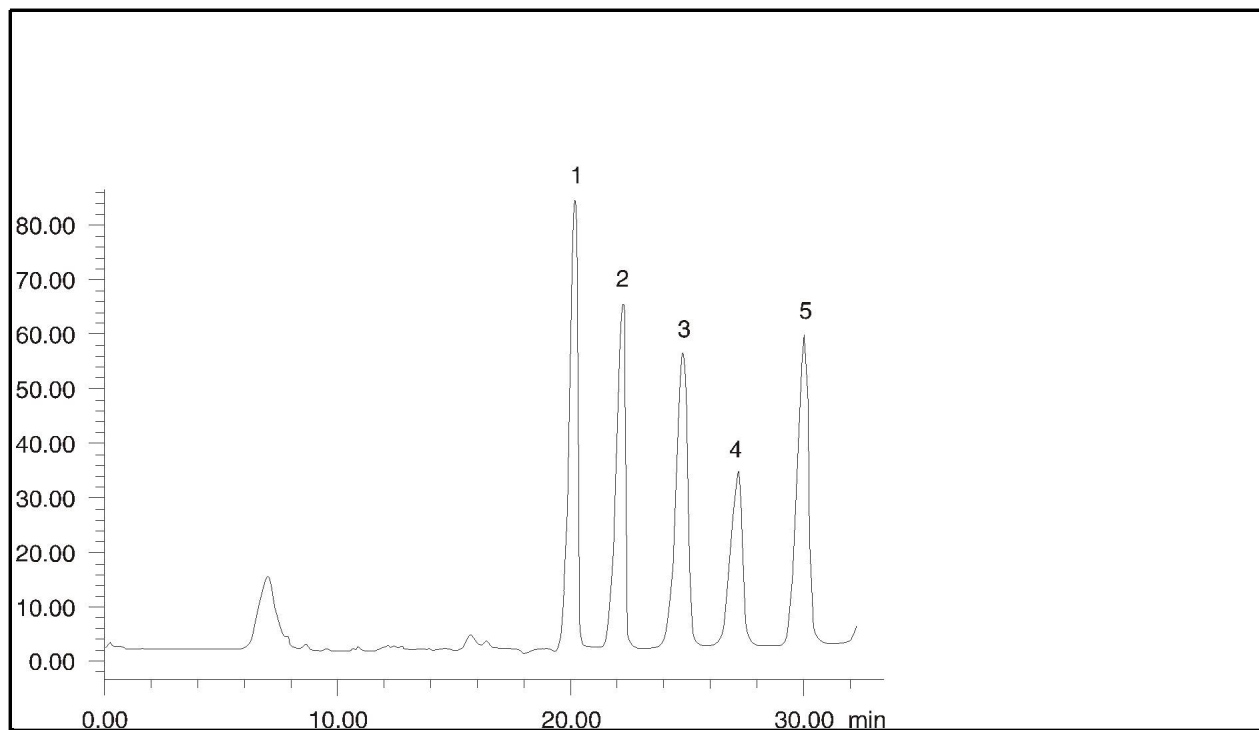

Column: Rezex 8 μ m Sulfonic Acid (SCX-SO₃H), strong cation exchange RCM-Monosaccharide

Dimensions: 300 x 7.8 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Water

Pct. A: 100

Flow Rate: 0.5 mL/min

Temperature: 60 °C

Detection: Refractive Index @ nm

SAMPLE ANALYTE

- 1 Glycerol
- 2 Methoxypolyethylene glycol
- 3 Triethylene glycol
- 4 Sorbitol
- 5 Urea

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 5518/ALF/537