

Carbamate/Urea Pesticides on Envirosep-PP (EPA Method 632)

Page 1 of 1



Column: EnviroSep-PP
Dimensions: 125 x 4.6 mm ID
Order No: 00E-3029-E0
Elution Type: Gradient
Eluent A: Water
Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	1.5	85	15
	3	11.5	0	100

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

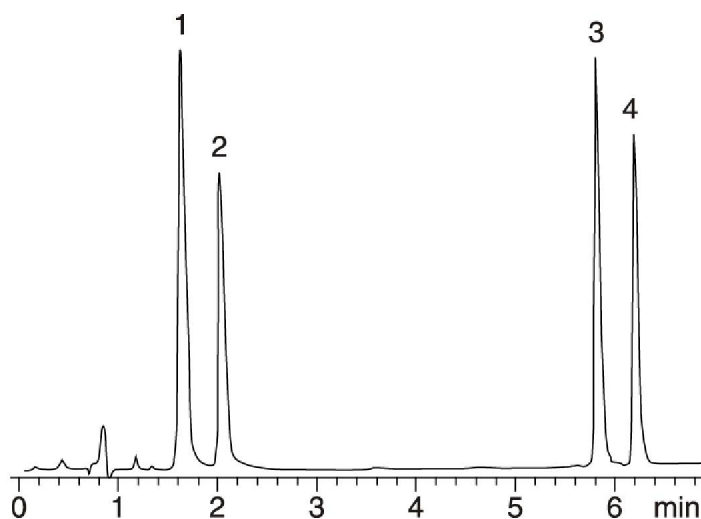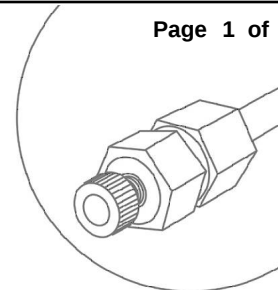

ANALYTES:

- 1 Oxamyl
- 2 Methomyl
- 3 Carbofuran
- 4 Fluometuron



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.