

Benzidines on Envirosep-PP (EPA Method 605)

Page 1 of 1

Column: EnviroSep-PP

Dimensions: 125 x 3.2 mm ID

Order No: 00E-3029-R0

Elution Type: Isocratic

Eluent A: 50 mM Sodium acetate/0.1% TEA, pH 4.7

Pct. A: 50



Eluent B: Acetonitrile

Pct. B: 50

Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

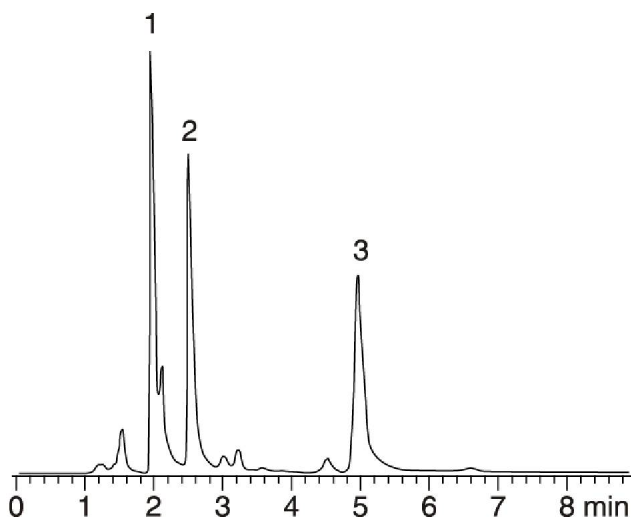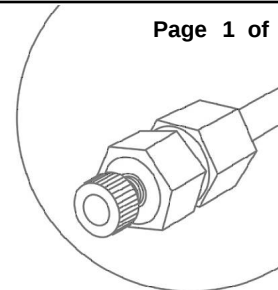

ANALYTES:

- 1 Benzidine
- 2 3,3'-Dimethylbenzidine
- 3 3,3'-Dichlorobenzidine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.