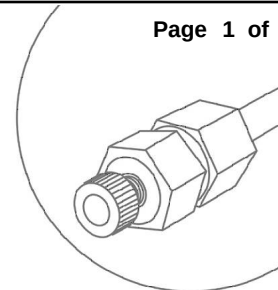
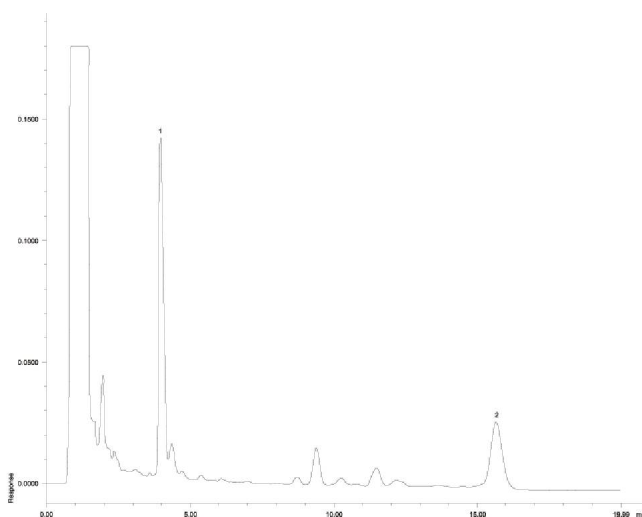


St. John's Wort (*Hypericum perforatum*)**Column:** Luna 5u C18(2) 100A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4252-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 50**Eluent B:** Methanol**Pct. B:** 30**Eluent C:** 0.1N Ammonium Acetate, pH 9.7**Pct. C:** 20**Flow Rate:** 1.5 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs. - Variable Wave. @ 330 nm (ambient)

Page 1 of 1

**Products used in this Application:****LUNA®****ANALYTES:**

- 1 Pseudohypericin
- 2 Hypericin

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.