

Fat-soluble Vitamins

Column: Luna 5u Silica 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4042-E0

Elution Type: Isocratic

Eluent A: n-heptane

Eluent B: 2-propanol


Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 294 nm (ambient)

Pct. A: 99.8


Pct. B: 0.2

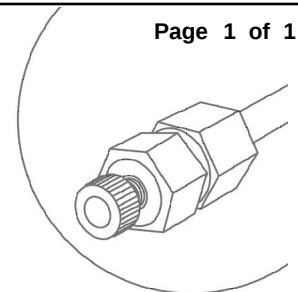

Page 1 of 1

Products used in this Application:

LUNA®

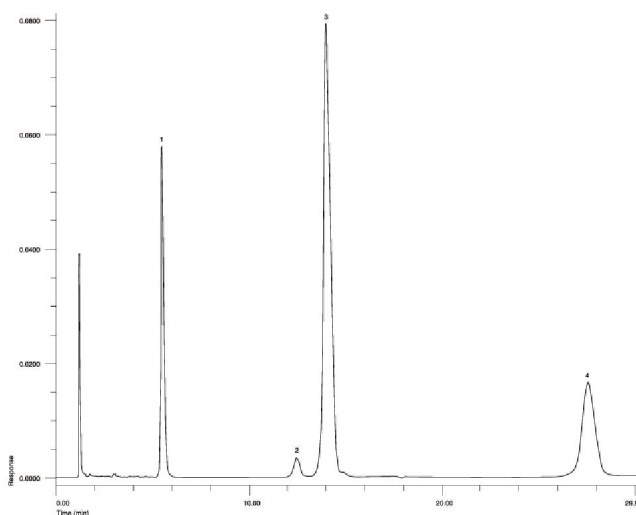

ANALYTES:

- 1 alpha-Tocopherol
- 2 beta-Tocopherol
- 3 gamma-Tocopherol
- 4 delta-Tocopherol



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.