



Proteins + Synthetic Peptides

Column: BioSep-SEC-S 2000**Dimensions:** 300 x 7.8 mm ID**Order No:** 00H-2145-K0**Elution Type:** Isocratic**Eluent A:** 50 mM Phosphate buffer, pH 6.8**Pct. A:** 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 215 nm (ambient)

Page 1 of 1

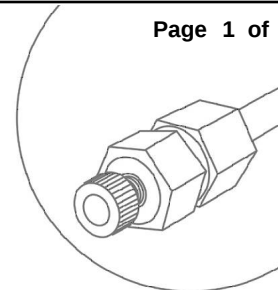

Products used in this Application:

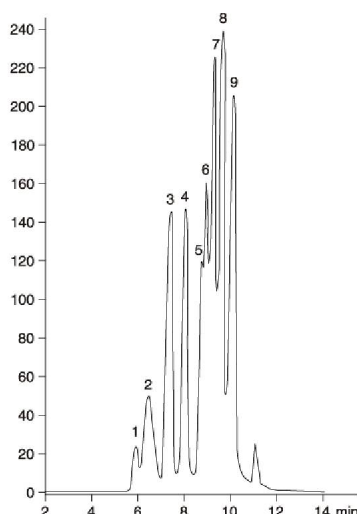Bio.Sep™

ANALYTES:

- 1 IgA (300,000)
- 2 IgG (150,000)
- 3 Ovalbumin (44,000)
- 4 Trypsin inhibitor (20,000)
- 5 Peptide + 5 (3,897)
- 6 Peptide + 4 (3,129)
- 7 Peptide + 3 (2,326)
- 8 Peptide + 2 (1,595)
- 9 Peptide + 1 (827)

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.