


Polymethyl Methacrylates (Wide MW Range)

Column: Phenogel 5u 10⁵Å
Dimensions: 300 x 7.8 mm ID
Order No: 00H-0446-K0
Elution Type: Isocratic
Eluent A: Tetrahydrofuran
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Differential refractometer (ambient)

Pct. A: 100


Page 1 of 1

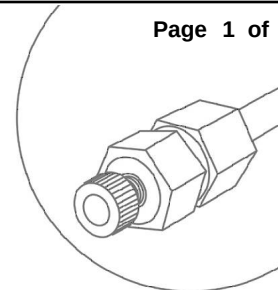

Products used in this Application:

Phenogel™

Analyst Note: 4 Phenogel 5μ columns (300 x 7.8mm ID) placed in series as follows: 10⁵Å, 10⁴Å, 10³Å, 500Å.

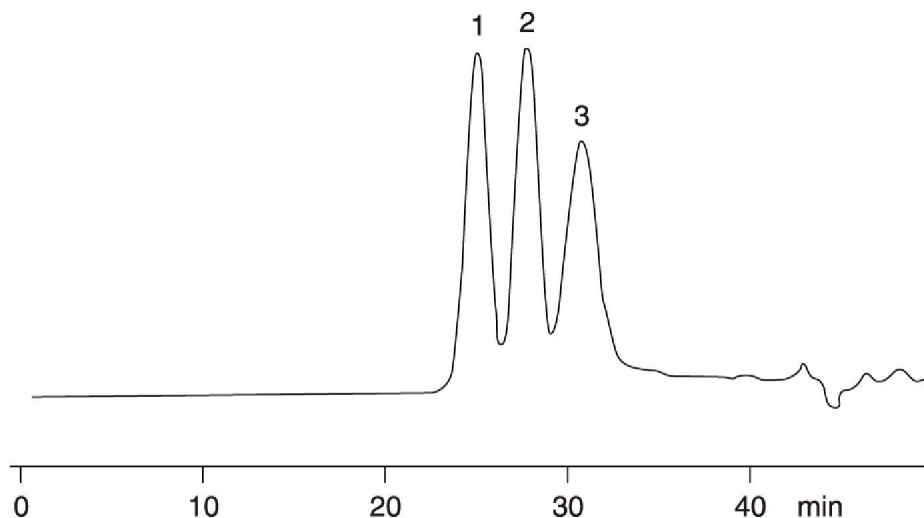

ANALYTES:

- 1 700,000 MW
- 2 107,000 MW
- 3 27,000 MW

