

Cyclic Polymer Characterization

Column: Phenogel 10u Linear(2)

Dimensions: 600 x 7.8 mm ID

Order No: 00K-3260-K0

Elution Type: Isocratic

Eluent A: Tetrahydrofuran with 1% TEA


Pct. A: 100

Flow Rate: 1.3 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 268 nm (ambient)

Page 1 of 1

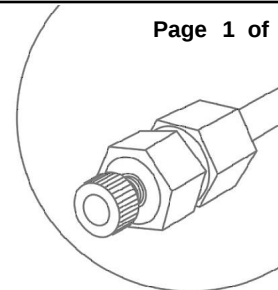

Products used in this Application:

Phenogel™

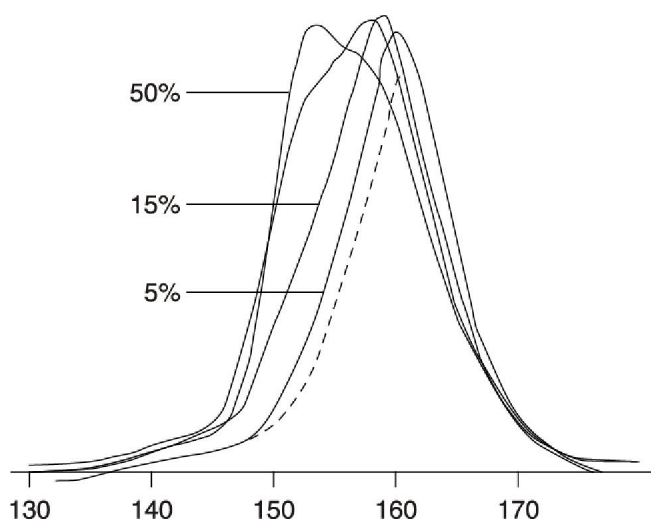

ANALYTES:

- 1 Poly-(2-vinylpyridine)



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.