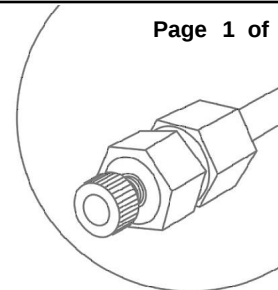
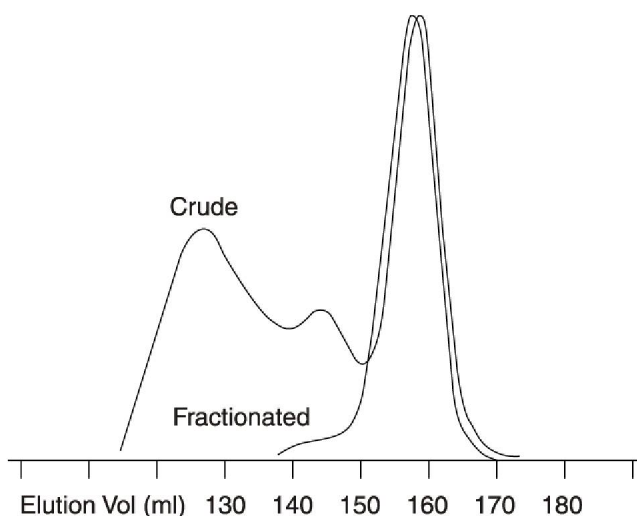


Cyclic Polymer Characterization

Column: Phenogel 10u Linear(2)**Dimensions:** 600 x 7.8 mm ID**Order No:** 00K-3260-K0**Elution Type:** Isocratic**Eluent A:** Tetrahydrofuran with 1% TEA**Pct. A:** 100**Flow Rate:** 1.3 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 268 nm (ambient)

Page 1 of 1

**Products used in this Application:****Phenogel[™]****ANALYTES:**

- 1 Poly-(2-vinylpyridine)

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.