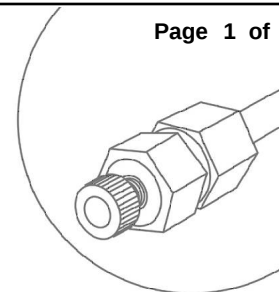
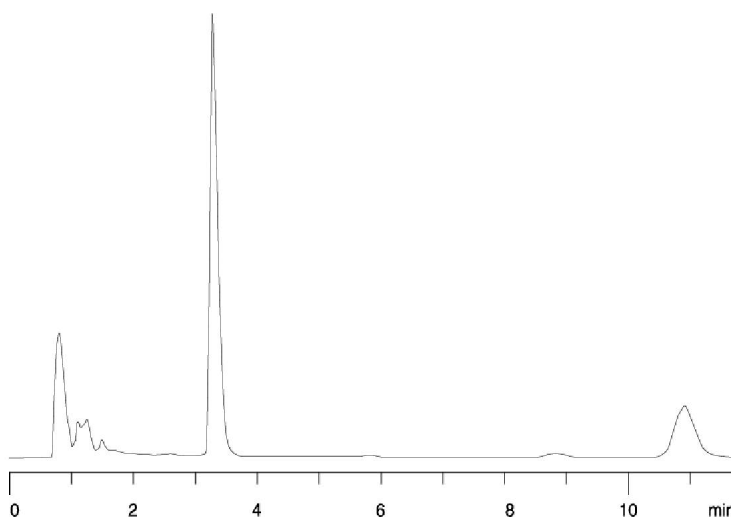


St. John's Wort on Synergi Hydro-RP

Page 1 of 1

Column: Synergi 4u Hydro-RP 80A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4375-E0**Elution Type:** Isocratic**Eluent A:** 100mM Ammonium acetate**Pct. A:** 20**Eluent B:** Acetonitrile**Pct. B:** 50**Eluent C:** Methanol**Pct. C:** 30**Flow Rate:** 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 590 nm (ambient)

Products used in this Application:

synergi[™]

ANALYTES:

- 1 Pseudohypericin
- 2 Hypericin



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.