


FAMES

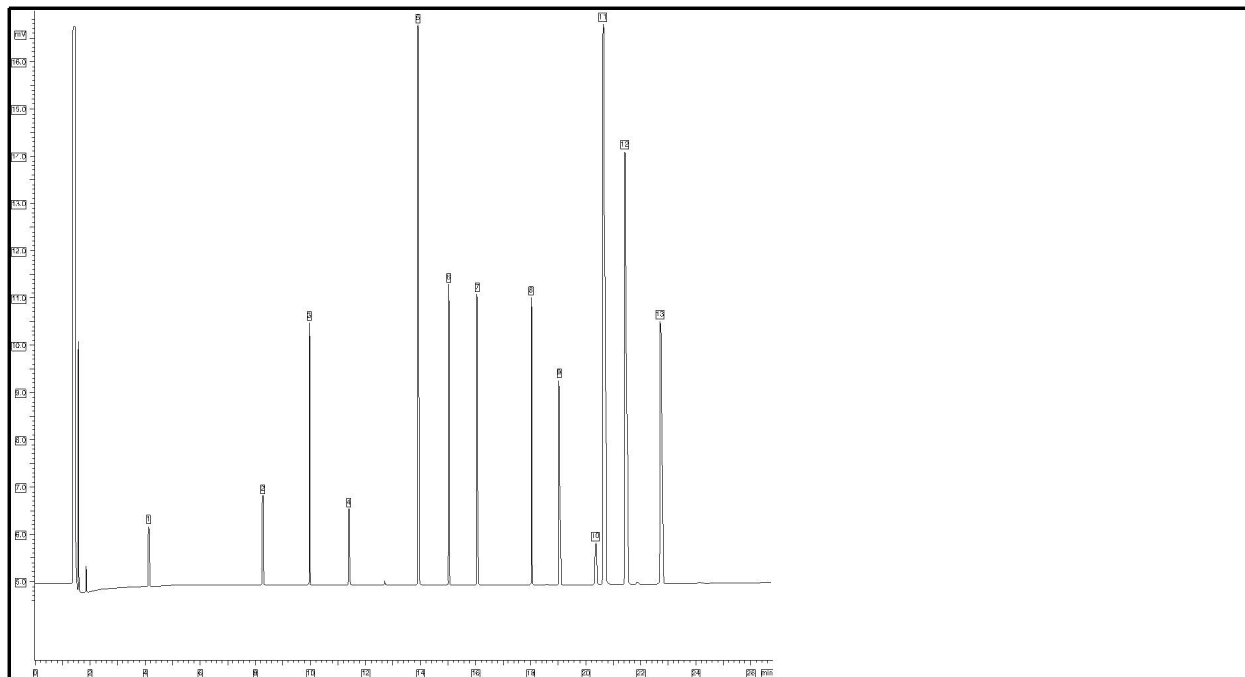

Column: Zebtron ZB-WAX
Phase: Polyethylene glycol (PEG)
Dimensions: 30 meters x 0.32 mm x 0.25 μ m
Order No:

Gradient Profile:

Step No.	Step Time	Flow Rate	Temperature
1	0	0	70
2	5	0	70
3	18	0	200
4	28	0	200

Carrier Gas: Helium 20.72 cm/sec
Flow Rate: 1 mL/min
Injection: Split 50:1 0 μ L @ 230 $^{\circ}$ C
Detection: Flame Ionization (FID) 250 $^{\circ}$ C

SAMPLE ANALYTE

1	Caproic acid	10	Stearic acid
2	Caprylic acid	11	Oleic acid
3	Nonanoic acid	12	Linoleic acid
4	Capric acid	13	Linolenic acid
5	Lauric acid		
6	Tridecanoic acid		
7	Myristic acid		
8	Palmitic acid		
9	Heptadecanoic acid		

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 6569/ZEB/016