

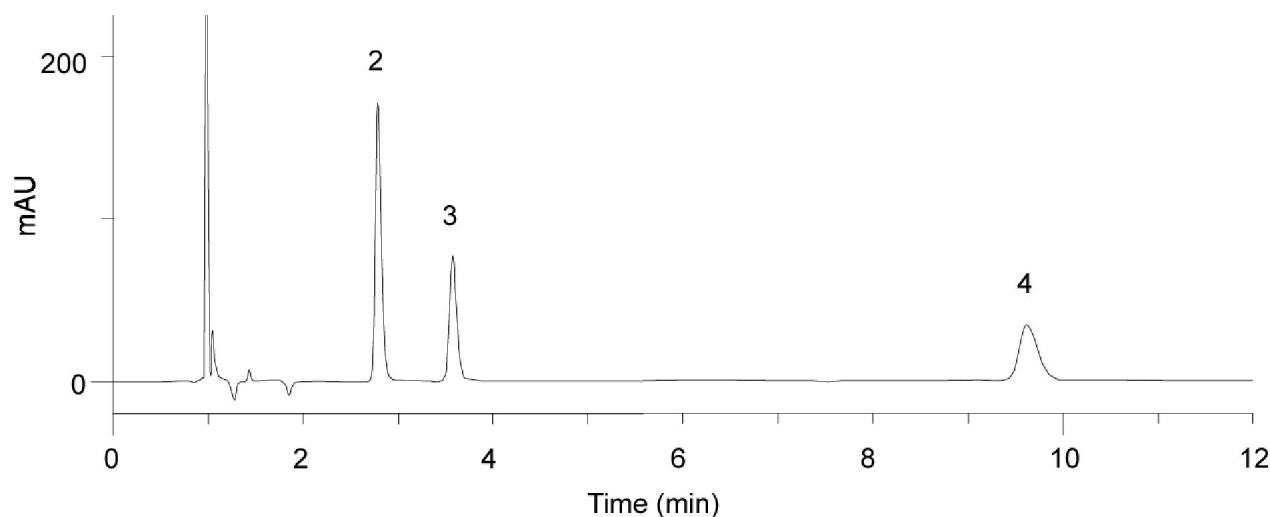
Cough and Colds on Synergi Hydro-RP

Page 1 of 1

Column: Synergi 4u Hydro-RP 80A
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4375-E0
Elution Type: Isocratic
Eluent A: Postassium Phosphate Buffer pH 2.5
Eluent B: Acetonitrile
Flow Rate: 1.5 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Pct. A: 75**Pct. B:** 25

Products used in this Application:

synergi™

ANALYTES:

- 1 Maleic Acid
- 2 Triprolidine
- 3 Chlorpheniramine
- 4 Diphenhydramine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.