

Cough and Cold on Luna SCX

Column: Luna 5u SCX 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4398-E0

Elution Type: Isocratic

Eluent A: 50mM Sodium Phosphate, pH 2.7

Pct. A: 40


Eluent B: Acetonitrile

Pct. B: 60

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

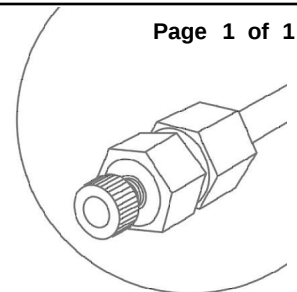

Page 1 of 1

Products used in this Application:

LUNA®

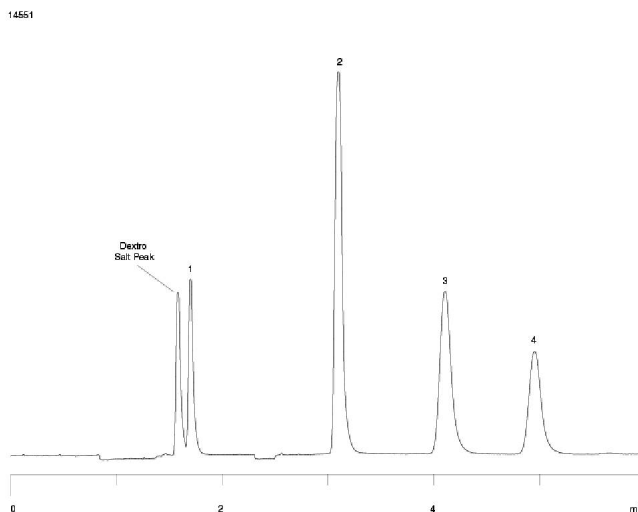

ANALYTES:

- 1 Benzoic Acid
- 2 Diphenhydramine
- 3 Dextromethorphan
- 4 Phenylephrine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.