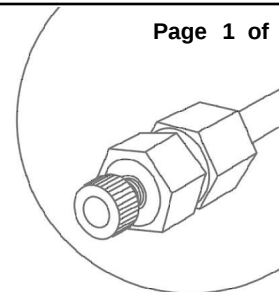
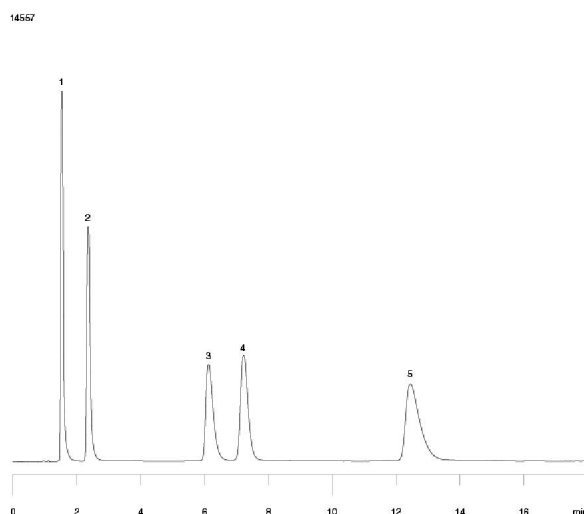


Deoxynucleosides on Luna SCX

Column: Luna 5u SCX 100A**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4398-E0**Elution Type:** Isocratic**Eluent A:** 50mM Potassium Phosphate, pH 2.5**Pct. A:** 100**Flow Rate:** 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Page 1 of 1

**Products used in this Application:****LUNA[®]****ANALYTES:**

- 1 d-Uridine
- 2 d-Inosine
- 3 d-Guanosine
- 4 d-Cytidine
- 5 d-Adenosine

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.