


Chiral Test Mix No. 4 (AL0-3049)

Column: Chirex 3126 (D)-penicillamine
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3126-E0
Elution Type: Isocratic
Eluent A: 2mM copper sulfate pentahydrate in H₂O/IPA (95:5) **Pct. A:** 100
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)


Page 1 of 1

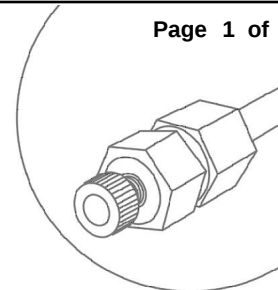

Products used in this Application:

CHIREX[®]

Analyst Note: This standard is used to test the following Chirex Chiral Columns
Phase: 3126.
NOTE: Separation conditions and results may vary for your column. Contact your Technical Representative.

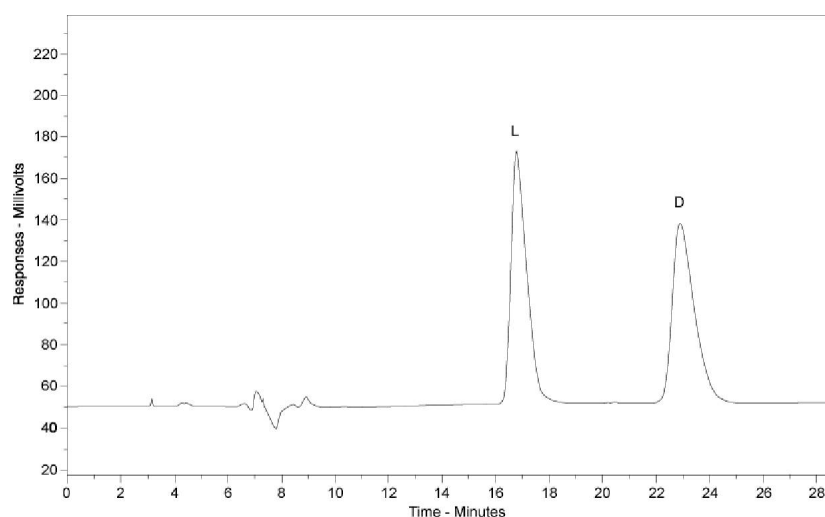

ANALYTES:

- 1 L-aspartic acid
- 2 D-aspartic acid

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.