

Proteins on Jupiter C4 300A

Column: Jupiter 5u C4 300A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4167-E0

Elution Type: Gradient

Eluent A: 0.1% TFA : H₂O

Eluent B: 0.08% TFA : ACN


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.5	95	5
	2	20	20	80

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. (ambient)

Page 1 of 1

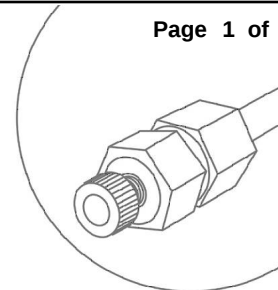

Products used in this Application:

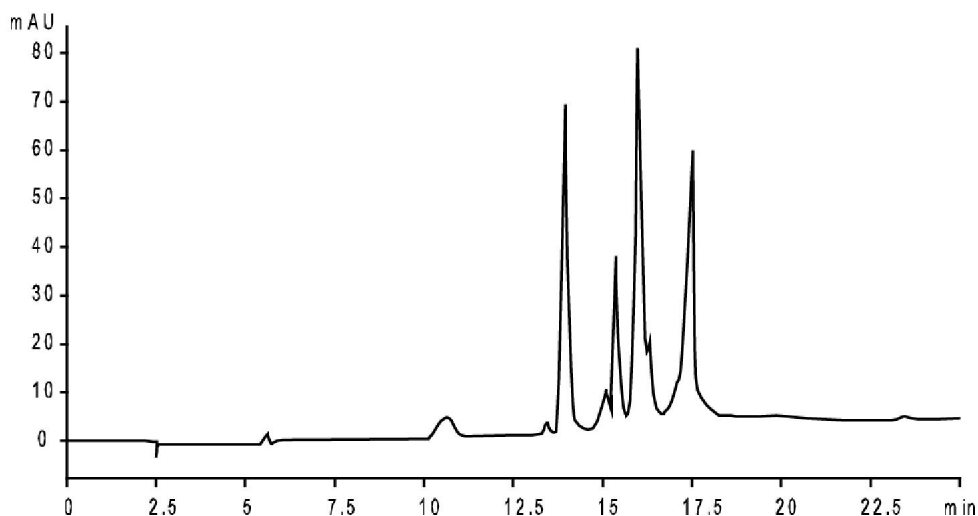

ANALYTES:

- 1 BSA
- 2 Glutamic Dehydrogenase
- 3 B-Galactosidase
- 4 Ovalbumin

