



EPA Method 8015: Nonhalogenated Organics

Page 1 of 1

Column: Zebron ZB-5
Phase: 5%-diphenyl-95%-dimethyl polysiloxane copolymer
Dimensions: 30 meters x 0.53 mm x 1.5 µm
Order No: 7HK-G002-28
Oven Profile: 45°C for 3 min to 60°C at 5°C/min to 275°C at 20°C/min
Carrier Gas: Constant Flow Helium, 4.8 mL/min
Injection: split 20:1 1 µL @ 225 °C
Detection: Flame Ionization (FID) (350 °C)

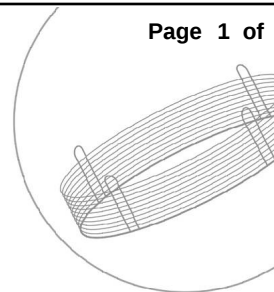

Products used in this Application:

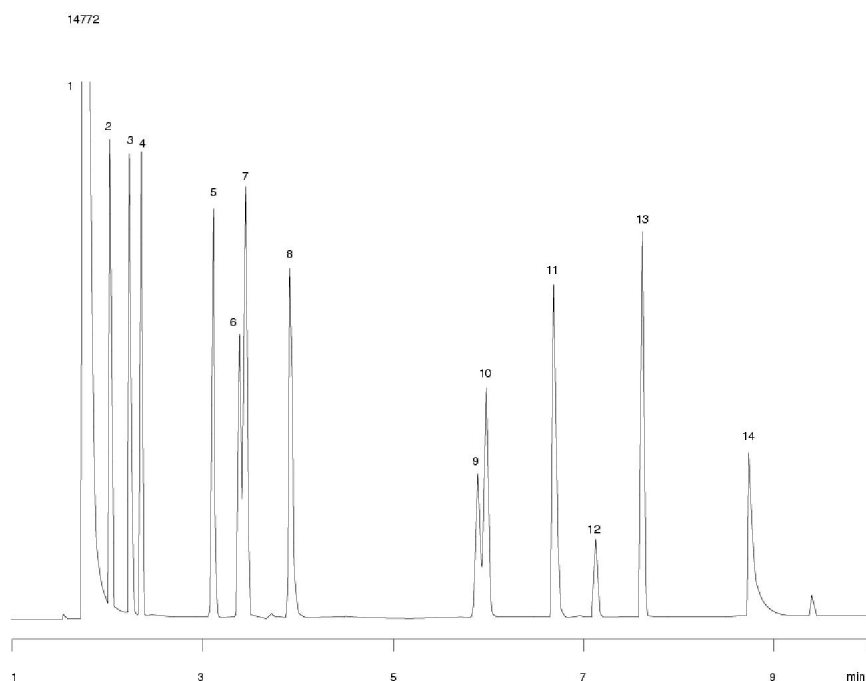Zebron™

ANALYTES:

- | | |
|---------------------------|--------------------------------|
| 1 Methanol (solvent) | 10 Methylmethacrylate |
| 2 Ethanol | 11 4-Methyl-2-Pentanone [MIBK] |
| 3 Acetonitrile | 12 Impurity |
| 4 Diethyl ether | 13 Ethylmethacrylate |
| 5 Propanenitrile | 14 Acrylamide |
| 6 2-Butanone [MEK] | |
| 7 2-methyl-2-propanitrile | |
| 8 Isobutanol | |
| 9 1,4-Dioxane | |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with tradition™

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au