


Aromatic Amines

Column: Zebron ZB-5MS
Phase: 5% Polysilarylene - 95% Polydimethylsiloxane copolymer
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 60°C for 1 min to 130°C at 12°C/min to 220°C at 30°C/min for 4 min
Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: split 100:1 1 µL @ 220 °C
Detection: Mass Selective (MSD) (240 °C)


Page 1 of 1

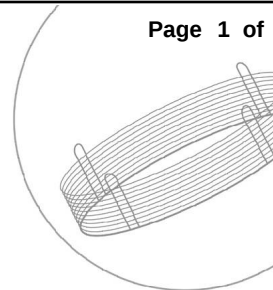

Products used in this Application:

Zebron™

15429

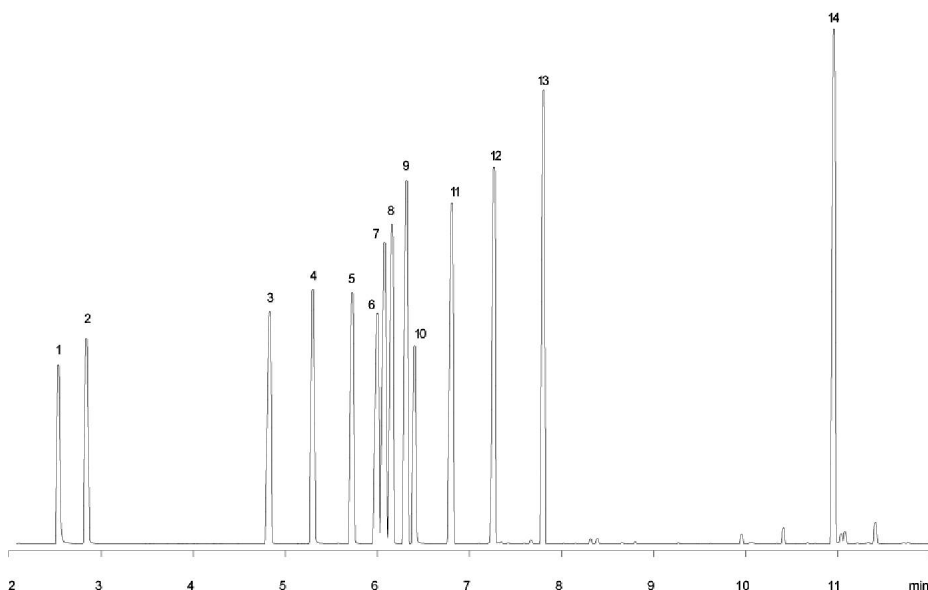

ANALYTES:

- | | |
|--------------------------|--------------------------|
| 1 Piperidine | 10 beta-Phenylethylamine |
| 2 2-Methylpiperidine | 11 N-Ethylaniline |
| 3 Aniline | 12 2,4-Dimethylaniline |
| 4 Benzylamine | 13 N,N-Diethylaniline |
| 5 alpha-Phenylethylamine | 14 Dibenzylamine |
| 6 N-methylaniline | |
| 7 m-Toluidine | |
| 8 o-Toluidine | |
| 9 N,N-dimethylaniline | |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with tradition™

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au