

Phospholipids in HILIC mode using Luna NH2

Page 1 of 1

Column: Luna 5u NH2 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4378-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	5	95
	2	15	50	50

Flow Rate: 1 mL/min


Col. Temp.: 30 °C

Detection: ELSD (ambient)

Products used in this Application:

LUNA[®]

15733

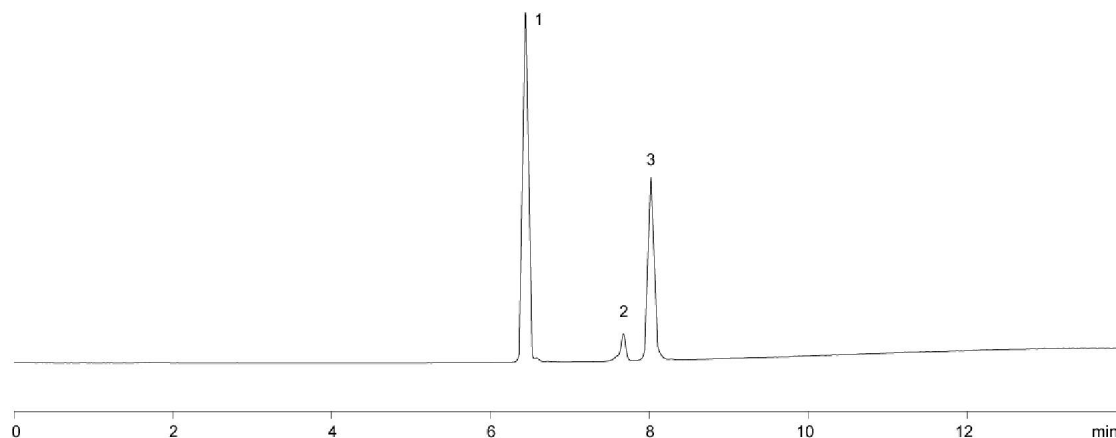

ANALYTES:

- 1 L- α -Phosphatidylcholine
- 2 L- α -Phosphatidylethanolamine
- 3 L- α -Phosphatidylinositol



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.