

## Basic Drugs at High pH on Gemini 5um C18

**Column:** Gemini 5u C18 110A

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4435-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM TEA pH11.0

**Pct. A:** 50

**Eluent B:** Acetonitrile


**Pct. B:** 50

**Flow Rate:** 1.5 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs. - Variable Wave. @ 230 nm (ambient)

Page 1 of 1

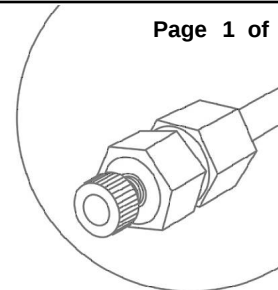

Products used in this Application:

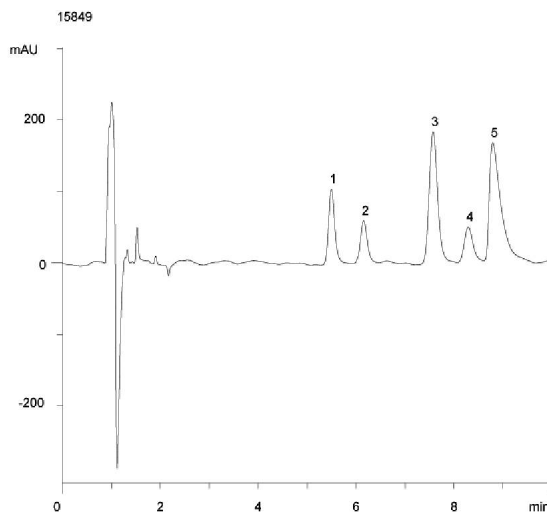

### ANALYTES:

- 1 Lidocaine
- 2 Tetracaine
- 3 Triprolidine
- 4 Protriptyline
- 5 Diphenhydramine

