


Non-Polar Basic Drugs by Gemini 5u C6-Phenyl

Page 1 of 1

Column: Gemini 5u C6-Phenyl 110A
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4444-E0
Elution Type: Isocratic
Eluent A: 20mM Potassium phosphate buffer, pH 2.5
Eluent B: Methanol
Flow Rate: 1.5 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Pct. A: 55

Pct. B: 45

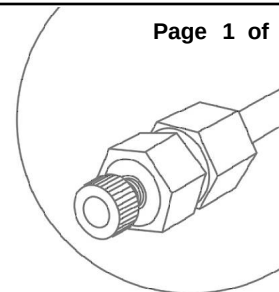


Products used in this Application:



Analyst Note: N (diphenhydramine) = 5533
 SYMM (diphenhydramine) = 0.65
 Backpressure = 217 bar

15997

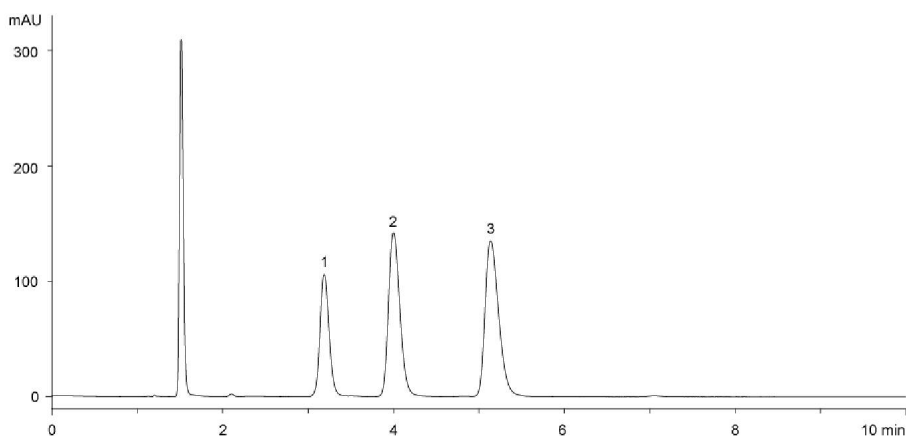

ANALYTES:

- 1 Chlorpheniramine
- 2 Triprolidine
- 3 Diphenhydramine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.