


Proteins

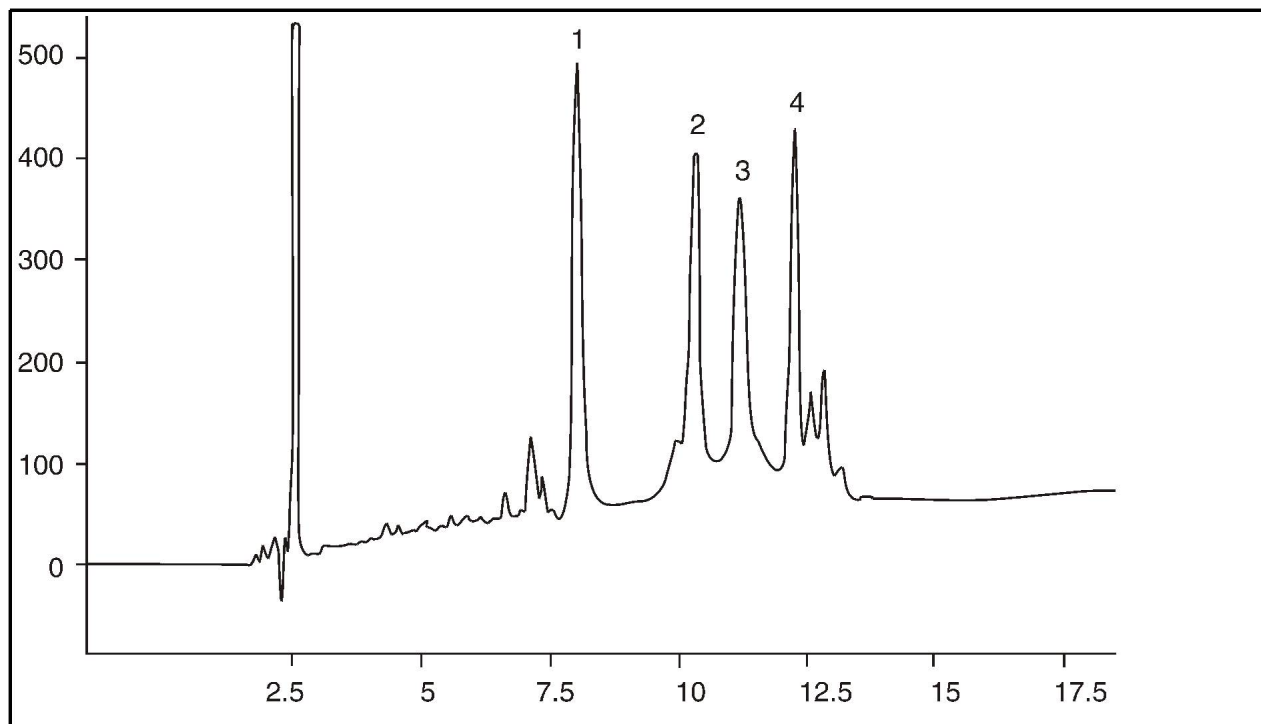

Column: Jupiter 5 μ m C5 300Å
Dimensions: 250 x 4.6 mm ID
Order No:
Elution Type: Gradient
Eluent A: 0.1% TFA in water
Eluent B: 0.1% TFA in acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	90	10
2	5	60	40

Flow Rate: 1 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ 220 nm

SPE:**SAMPLE ANALYTE**

- 1 *Bovine serum albumin*
- 2 *Peroxidase (horseradish)*
- 3 *Catalase (beef liver)*
- 4 *Carbonic anhydrase (bovine)*

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 8722/BRO/13