


Saponins

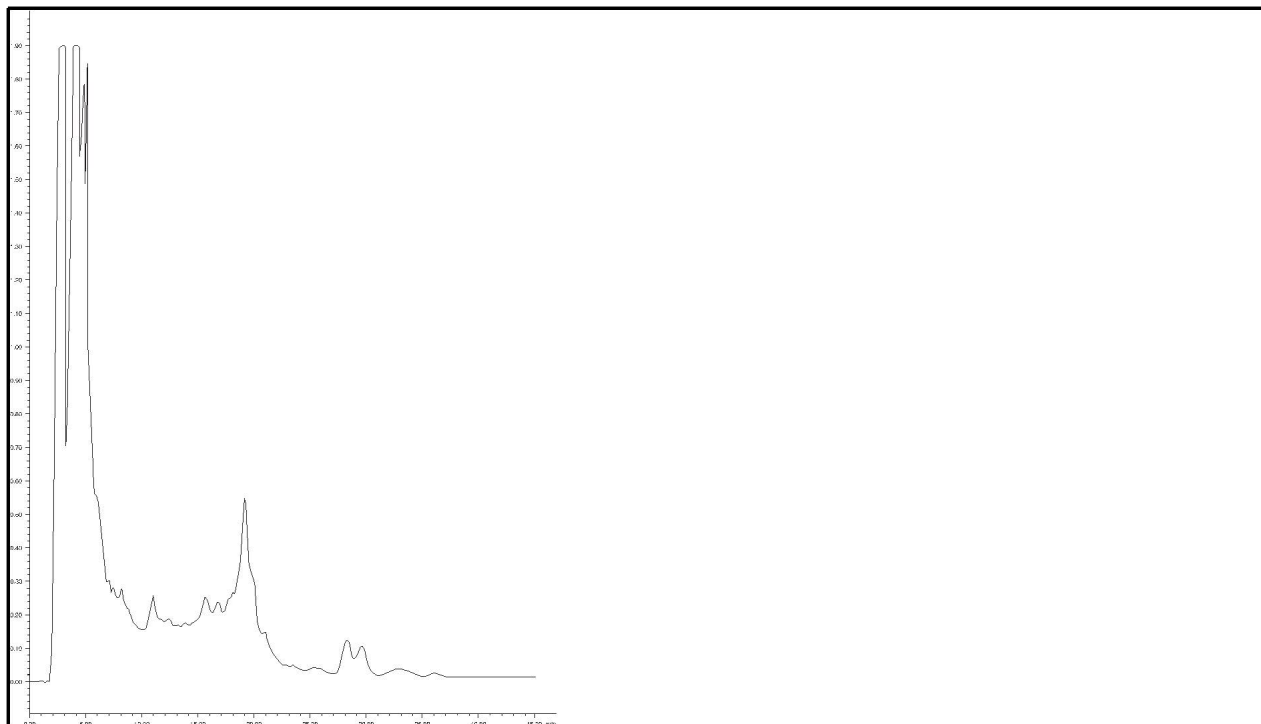

Column: Prodigy 10 μ m C18 / ODS Prep

Dimensions: 250 x 10 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile (various; 36% - 43%)

Pct. A: 43

Eluent B: Water

Pct. B: 57

Flow Rate: 4 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 203 nm

SAMPLE ANALYTE

1 Saponin

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 10128/PRO/036