


Analgesics and Anti-Inflammatories

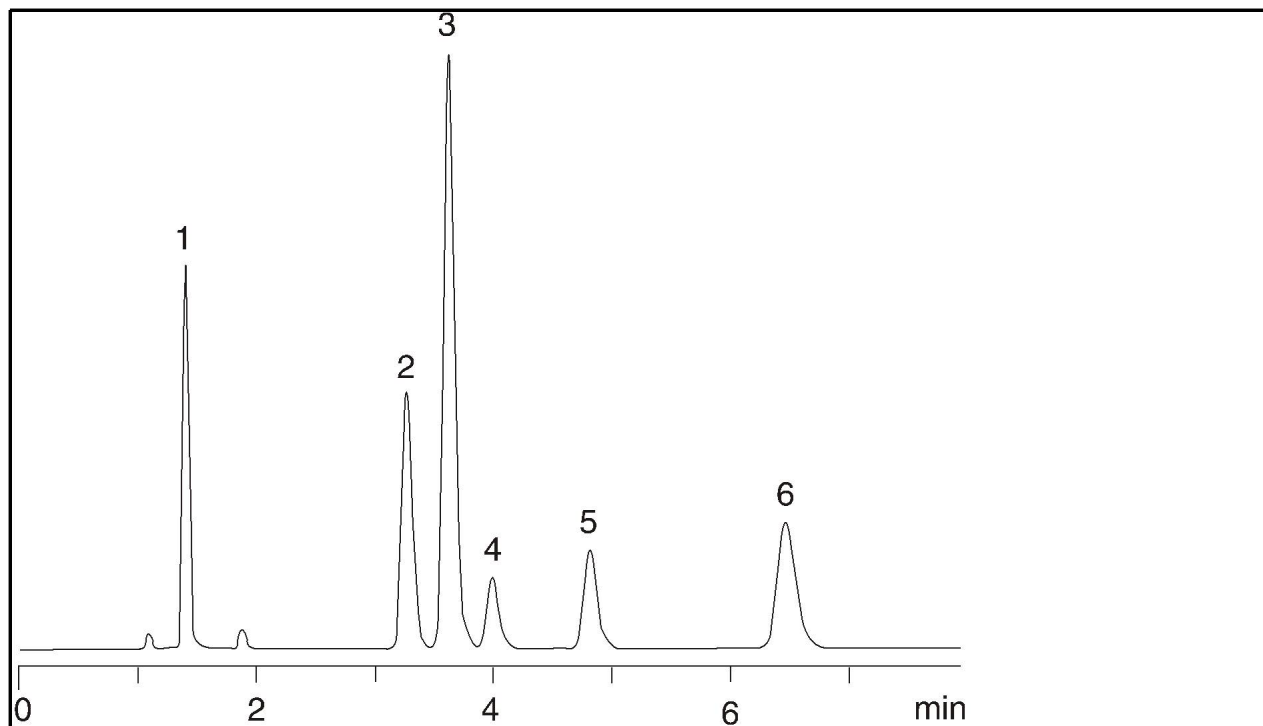

Column: Prodigy 5 μ m Phenyl PH-3

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20 mM Potassium phosphate pH 2.5

Pct. A: 50

Eluent B: Acetonitrile

Pct. B: 40

Eluent C: Methanol

Pct. C: 10

Flow Rate: 1.5 mL/min

Temperature: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 230 nm

SAMPLE ANALYTE

- 1 Acetaminophen
- 2 Indoprofen
- 3 Naproxen
- 4 Diflunisal
- 5 Ibuprofen
- 6 Indomethacin

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. 10254/BRO/5
No.