

Nucleotides

Column: Synergi 4u Polar-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate pH 2.7


Pct. A: 100

Flow Rate: 2 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

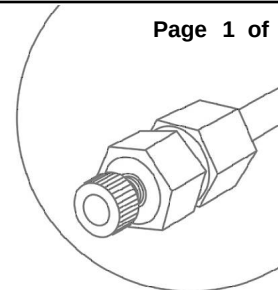

Products used in this Application:

synergi™

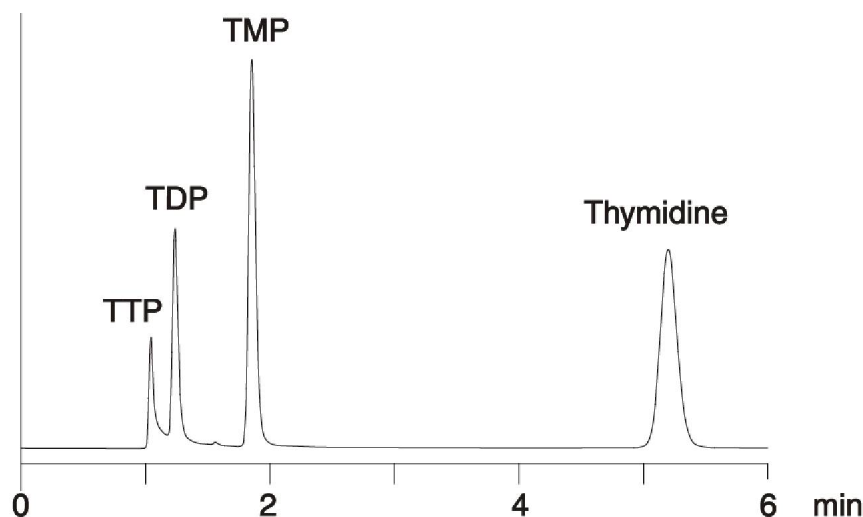

ANALYTES:

- 1 Thymidine triphosphate
- 2 Thymidine diphosphate
- 3 Thymidine monophosphate
- 4 Thymidine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.