


Opiates (BSTFA Derivatization)

Column: ZB-5
Dimensions: 30m x 0.25mm x 0.25µm
Flow: Constant flow of 1.1 mL/min. of He
Injector: 250°C; 1 µL at 10:1 split
Oven: 180°C to 230°C at 15°C/min. to 300°C at 4°C/min.
Detector: Scan mode MS
Analytes: All analytes were 6.6 ppm before derivatization except ethylmorphine at 666 ppb.

Opiates on ZB-5, 30m x 0.25mm x 0.25µm

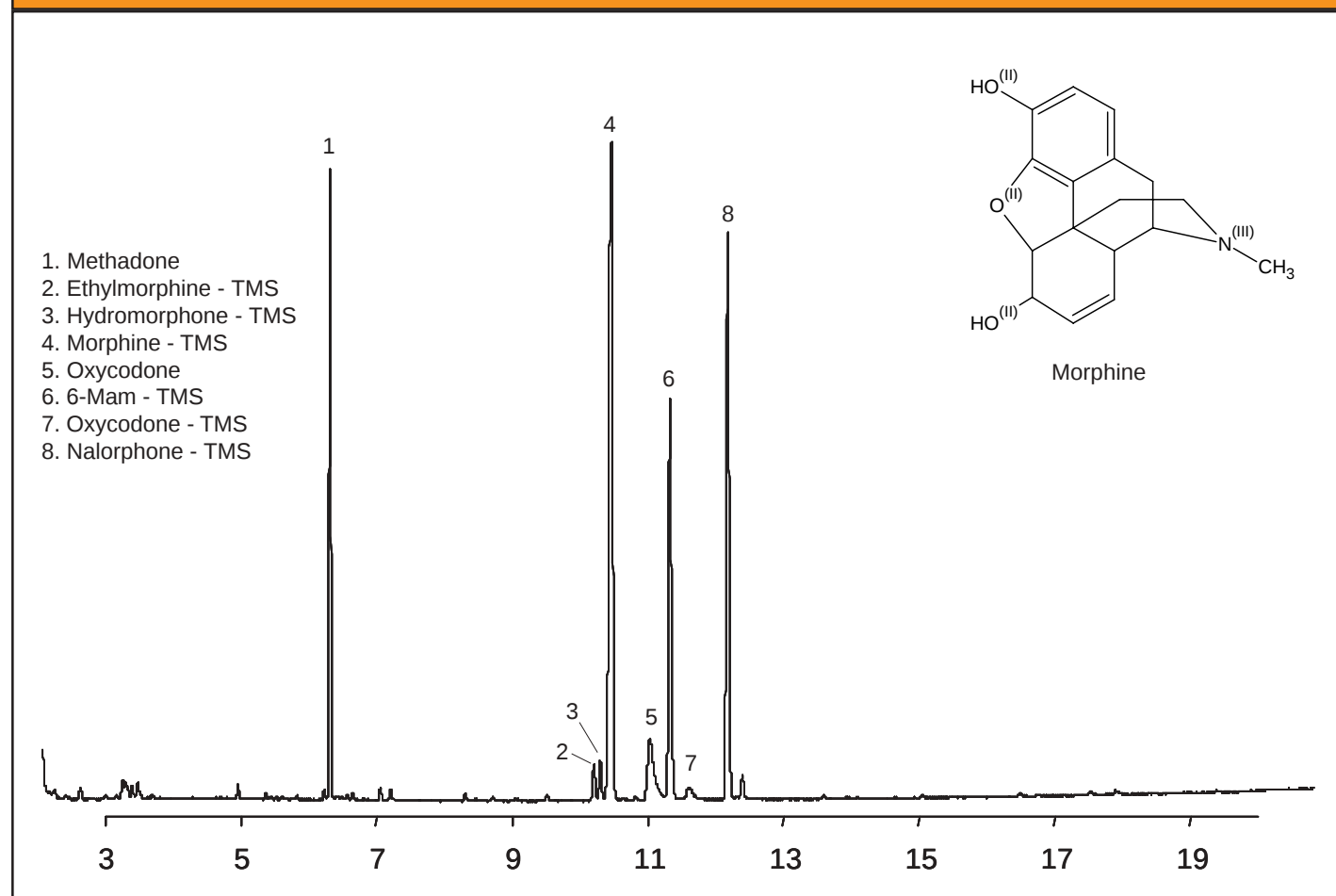

Ordering Information

Order Number	Description
7HG-G002-11	ZB-5 - 30m x 0.25mm x 0.25µm

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA	Puerto Rico	Canada	United Kingdom	Germany	New Zealand	Australia
tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	(800) 541-HPLC (310) 328-7768 info@phenomenex.com	(800) 543-3681 (310) 328-7768 info@phenomenex.com	01625-501367 01625-501796 ukinfo@phenomenex.com	06021-58830-0 06021-58830-11 anfrage@phenomenex.com	09-4780951 09-4780952 info@phenomenex.co.nz	1800-553-929 1800-553-923 info@phenomenex.com.au

phenomenex
...breaking with tradition