


Nucleotides

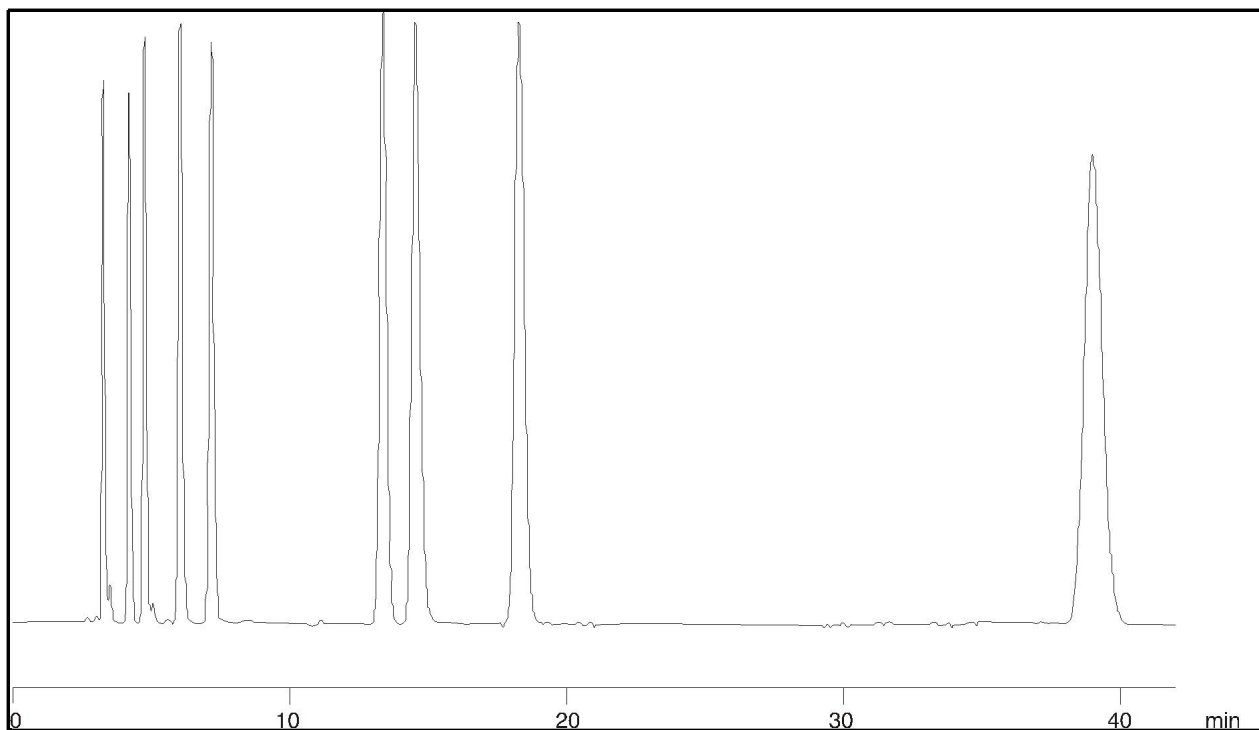

Column: Develosil 5 μ m C30 RP-Aqueous (C30-UG-5)

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 0.1 M KH₂PO₄ with KOH, pH 6.0

Pct. A: 100

Flow Rate: 1 mL/min

Temperature: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 260 nm

SAMPLE ANALYTE

- 1 5'-CTP
- 2 5'-CDP
- 3 5'-CMP
- 4 5'-GTP
- 5 5'-GDP
- 6 5'-GMP
- 7 5'-ATP
- 8 5'-ADP
- 9 5'-AMP

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 13551/A-4