


Tocopherols

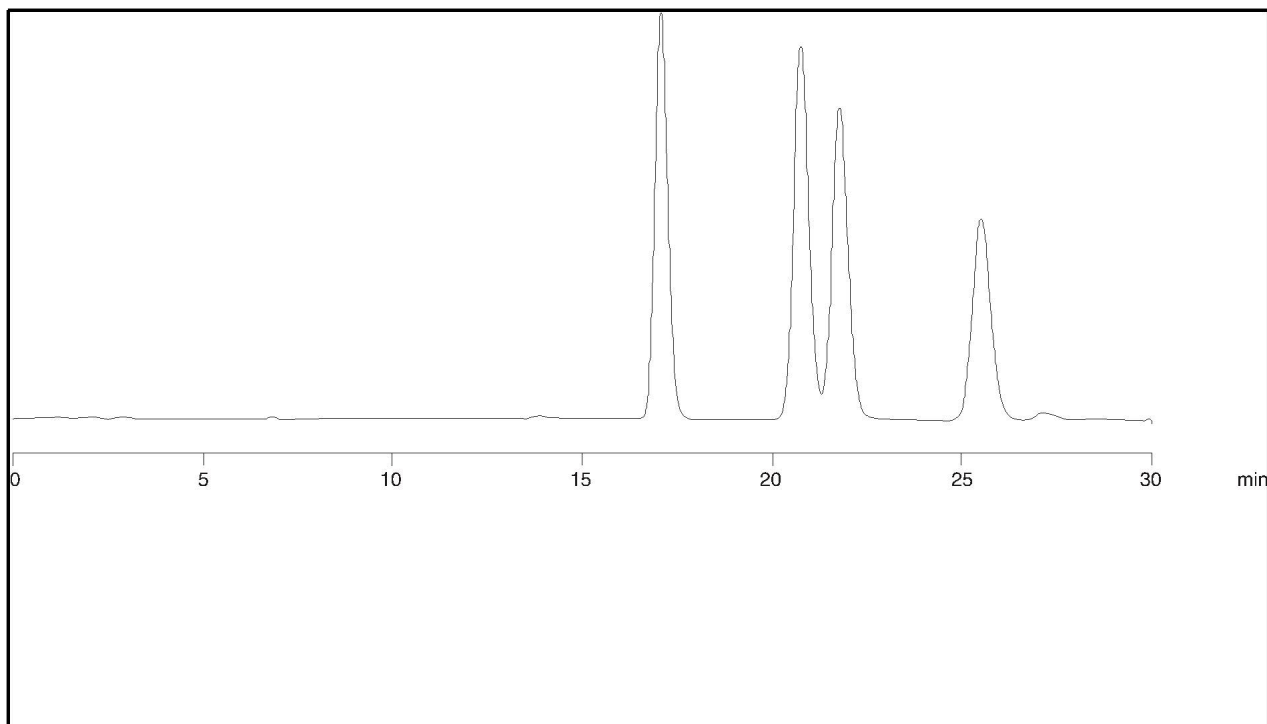

Column: Develosil 5 μ m C30 RP-Aqueous (C30-UG-5)

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Methanol/Water 97:3

Pct. A: 100

Flow Rate: 1 mL/min

Temperature: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 295 nm

SAMPLE ANALYTE

- 1 *delta-Tocopherol*
- 2 *gamma-Tocopherol*
- 3 *beta-Tocopherol*
- 4 *alpha-Tocopherol*

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 13574/A-14