

Antioxidants - Bondclone vs. μ Bondapak Comparison

Column: Bondclone 10u C18 148A

Dimensions: 300 x 3.9 mm ID

Order No: 00H-2117-C0

Elution Type: Gradient

Eluent A: Acetonitrile


Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	10	100	0


Flow Rate: 1.7 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

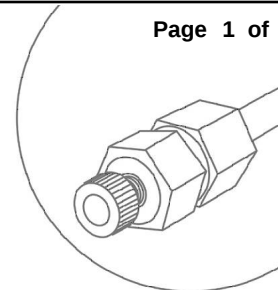

Products used in this Application:

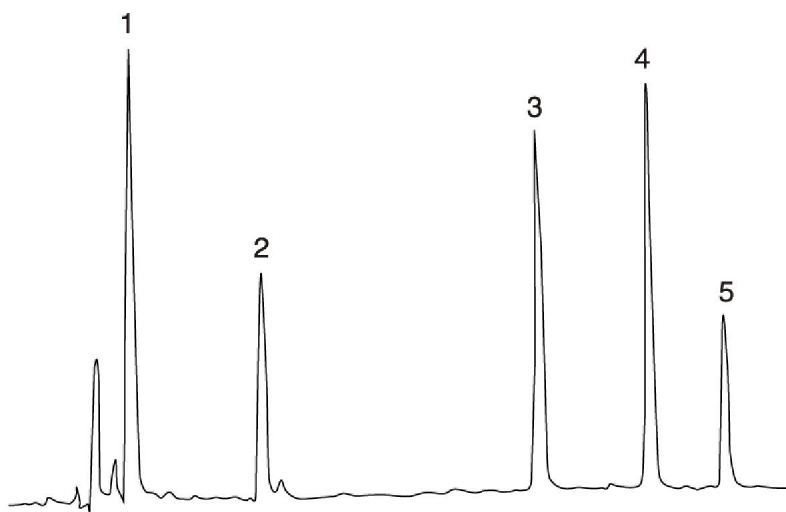Bondclone[™]

ANALYTES:

- 1 Ethyl 754
- 2 BHT
- 3 Goodrite 3125
- 4 Irganox 1010
- 5 Irganox 1076



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.