



## Uniconazole

**Column:** Chirex 3014 (S)-valine and (R)-1-( $\alpha$ -naphthyl)ethylamine**Dimensions:** 250 x 4 mm ID**Order No:** 00G-3014-D0**Elution Type:** Isocratic**Eluent A:** Hexane**Pct. A:** 82.6**Eluent B:** 1,2-Dichloroethane**Pct. B:** 16.5**Eluent C:** Ethanol**Pct. C:** 0.8**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

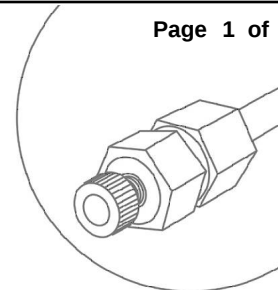

Products used in this Application:

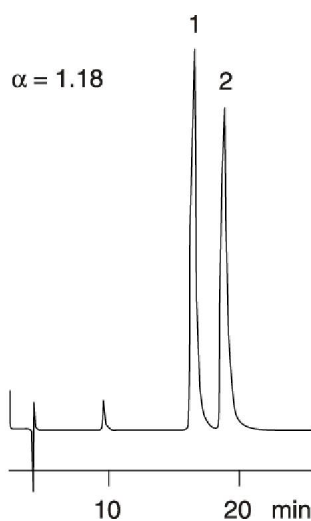**CHIREX<sup>®</sup>**

## ANALYTES:

- 1 R-Uniconazole
- 2 S-Uniconazole



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.