

DMSO and Ethanol

Column: Rezex ROA-Organic Acid H+ (8%)

Dimensions: 150 x 7.8 mm ID

Order No: 00F-0138-K0

Elution Type: Isocratic

Eluent A: 0.005 M H₂SO₄

Pct. A: 90

Eluent B: Acetonitrile


Pct. B: 10

Flow Rate: 0.6 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Page 1 of 1

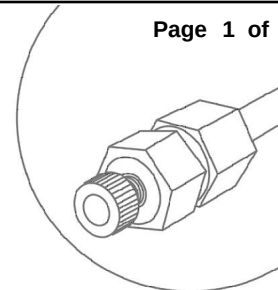

Products used in this Application:

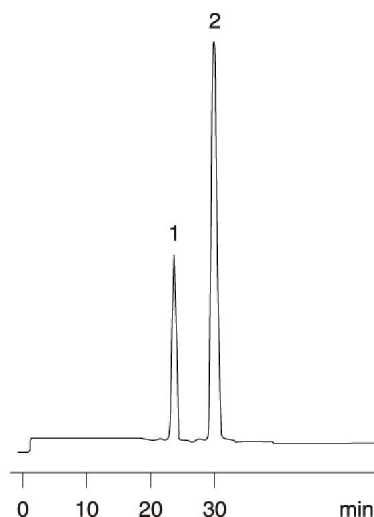REZEX[™]

ANALYTES:

- 1 Ethanol
- 2 DMSO

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.