

Ibuprofen

Column: CHIREX 3012 (R)-phenylglycine and 3,5-dinitroaniline

Dimensions: 250 x 4 mm ID

Order No: 00G-3012-D0

Elution Type: Isocratic

Eluent A: 0.01 M Ammonium acetate in methanol


Pct. A: 100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

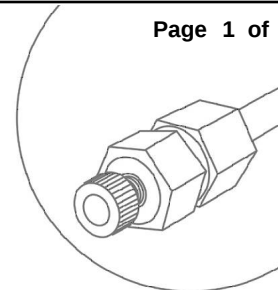

Products used in this Application:

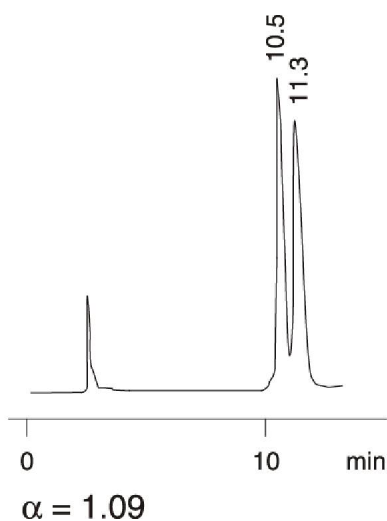CHIREX™

ANALYTES:

- 1 Ibuprofen
- 2 Ibuprofen

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.