

Diniconazole

Column: Chirex 3014 (S)-valine and (R)-1-(alpha-naphthyl)ethylamine

Dimensions: 250 x 4 mm ID

Order No: 00G-3014-D0

Elution Type: Isocratic

Eluent A: Hexane

Pct. A: 82.6

Eluent B: 1,2-Dichloroethane

Pct. B: 16.5

Eluent C: Ethanol


Pct. C: 0.8

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Page 1 of 1

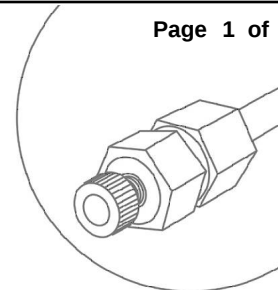

Products used in this Application:

CHIREX

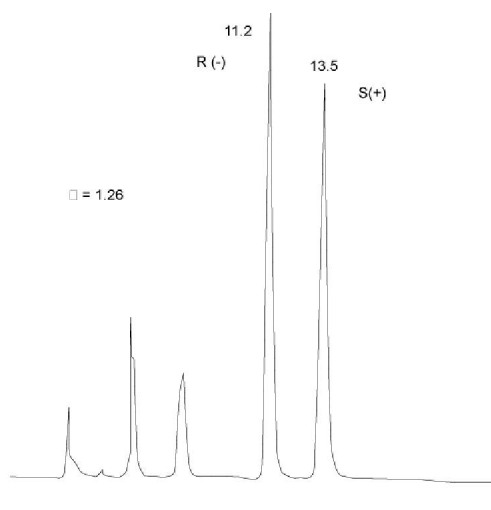

ANALYTES:

- 1 R(-)-Diniconazole
- 2 S(+)-Diniconazole



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.