



Derivatized amino acids in human urine

Page 1 of 1

Column: EZ:faast GC
Phase: Zebron ZB-AAA (proprietary)
Dimensions: 10 meters x 0.25 mm x 0 µm
Order No: KG0-7165
Oven Profile: 30°C/min from 110° to 320°C, hold at 320°C for 1 minute
Carrier Gas: Constant Flow Helium, 1.5 mL/min
Injection: Split 15:1 100 µL @ 250 °C
Detection: Flame Ionization (FID) (320 °C)

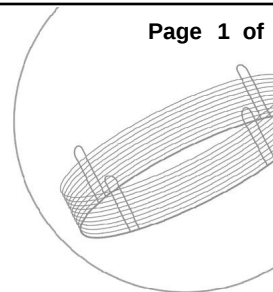

Products used in this Application:

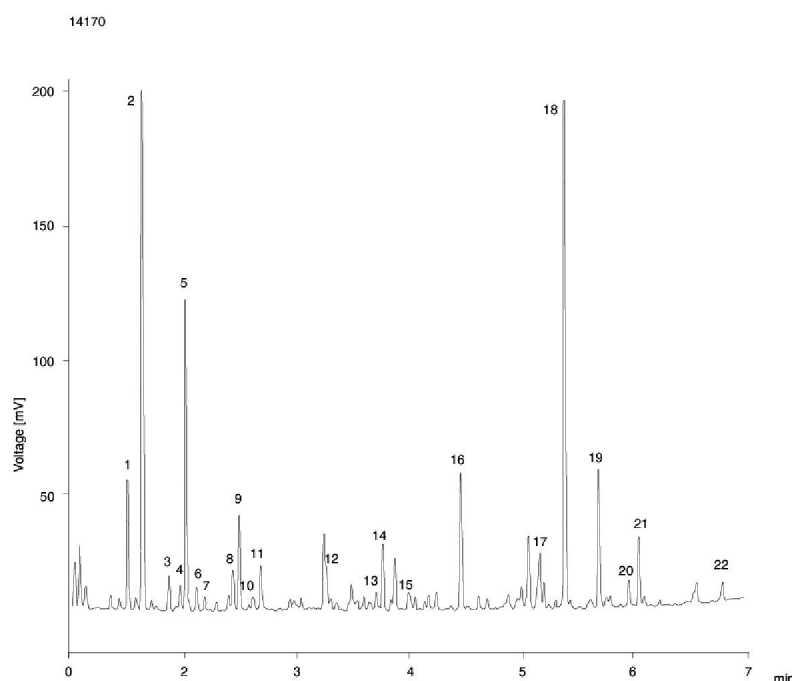EZ:faast[™]

ANALYTES:

- | | | |
|-----------------------------|---------------------------|---------------------------|
| 1 Alanine | 10 Proline | 19 Tyrosine |
| 2 Glycine | 11 Asparagine | 20 Prolinylhydroxyproline |
| 3 Valine | 12 Aspartic acid | |
| 4 Beta-aminoisobutyric acid | 13 Glutamic acid | |
| 5 Norvaline | 14 Phenylalanine | |
| 6 Leucine | 15 Alpha-aminoadipic acid | |
| 7 Isoleucine | 16 Glutamine | |
| 8 Threonine | 17 Lysine | |
| 9 Serine | 18 Histidine | |

©Copyright 2000-2005 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex[®]
...breaking with tradition[™]

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au