


American Ginseng Root (*Panax quinquefolium*)

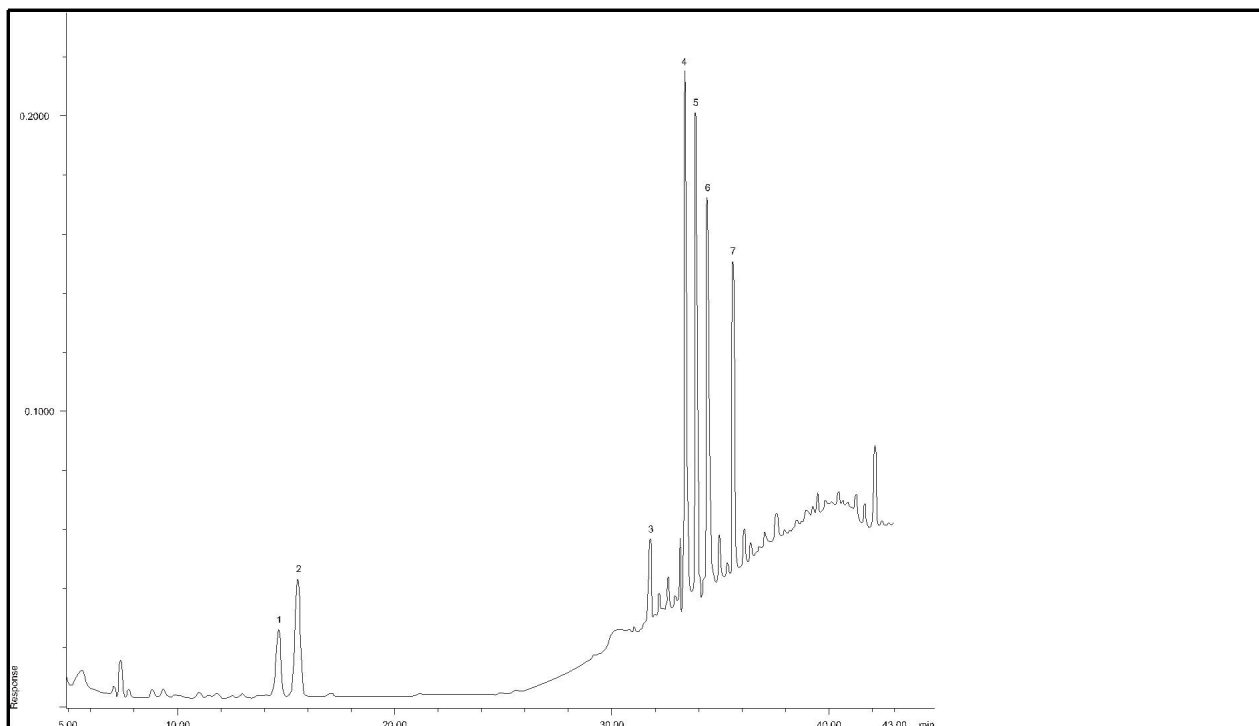

Column: Luna 3 μ m C18 / ODS (2)

Dimensions: 150 x 3 mm ID

Order No:

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	20	80
2	22	22	78
3	40	50	50
4	50	55	45
5	52	65	35
6	54	20	80
7	70	20	80

Flow Rate: 0.6 mL/min

Temperature: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 205 nm

SPE:

SAMPLE ANALYTE

- 1 Ginsenoside Rg 1
- 2 Ginsenoside Re
- 3 Ginsenoside Rf

- 4 Ginsenoside Rb 1
- 5 Ginsenoside Rc
- 6 Ginsenoside Rb 2


For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1123/NUT/ginsng2.doc

Continued...

American Ginseng Root (*Panax quinquefolium*)

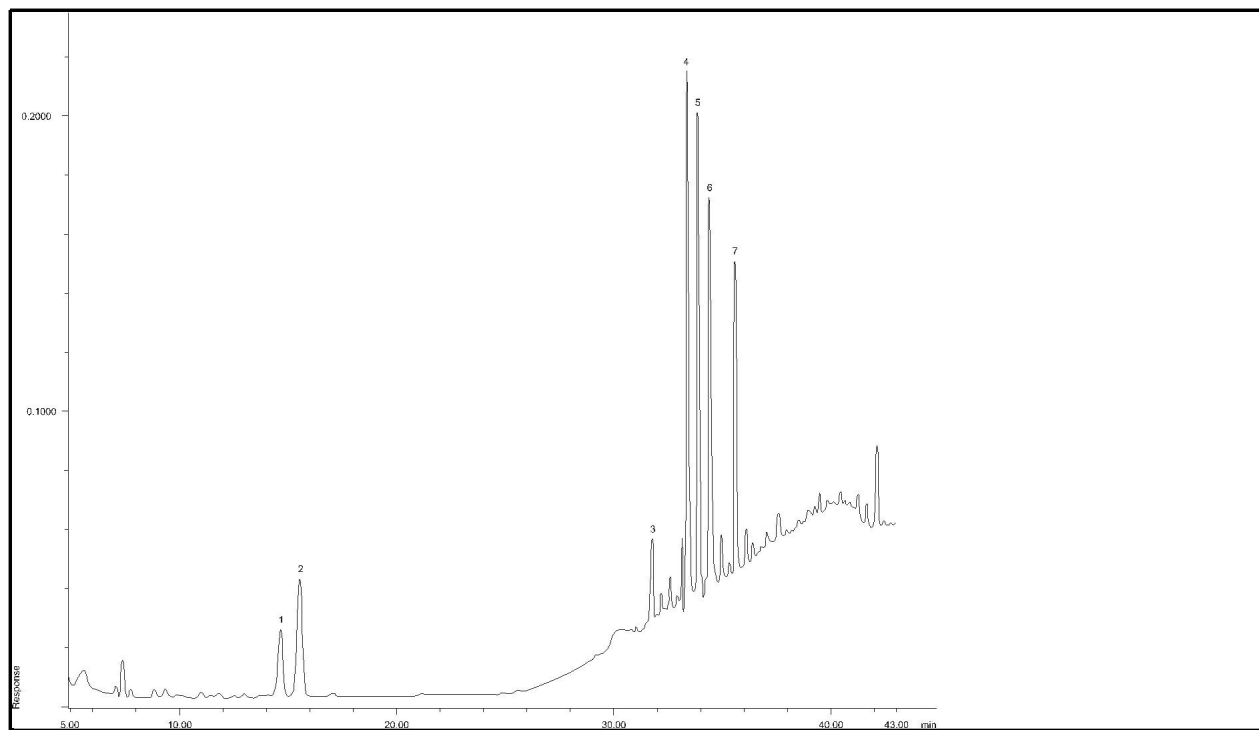

7 Ginsenoside Rd

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1123/NUT/ginsng2.doc