


Fat-soluble vitamins on Synergi Hydro-RP

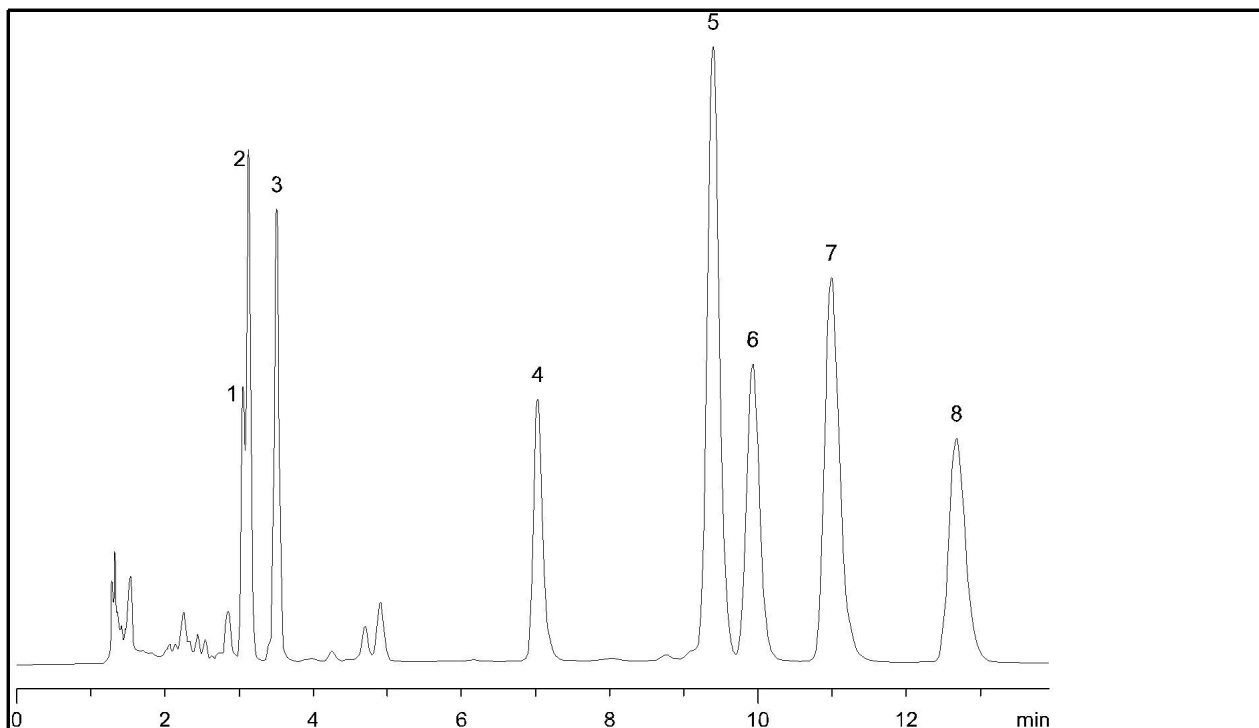

Column: Synergi 4 μ m C18 / ODS Hydro-RP

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 75

Eluent B: Methanol

Pct. B: 25

Flow Rate: 2 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm

SAMPLE ANALYTE

- 1 Impurity in Trans-Retinol (Vitamin A)
- 2 Vitamin A (trans-Retinol)
- 3 Vitamin A acetate
- 4 Vitamin K2
- 5 Vitamin D2 (Ergocalciferol)
- 6 Vitamin D3 (Cholecalciferol)
- 7 Vitamin E (alpha Tocopherol)
- 8 Vitamin E acetate (Tocopherol acetate)

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14214