



EPA Method 8081 (Organochlorine Pesticides) on ZB-1701P

Page 1 of 1

Column: Zebron ZB-1701P
Phase: 14%-cyanopropylphenyl-86%-dimethyl polysiloxane copolymer
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G012-11
Oven Profile: 100°C to 200°C at 25°C/min to 240°C at 6°C/min to 285°C at 20°C/min for 5 min
Carrier Gas: Constant Flow Helium, 1.6 mL/min
Injection: Splitless 1 µL @ 240 °C
Detection: Electron Capture (ECD) (300 °C)

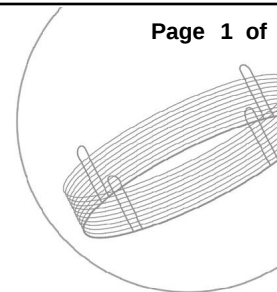

Products used in this Application:

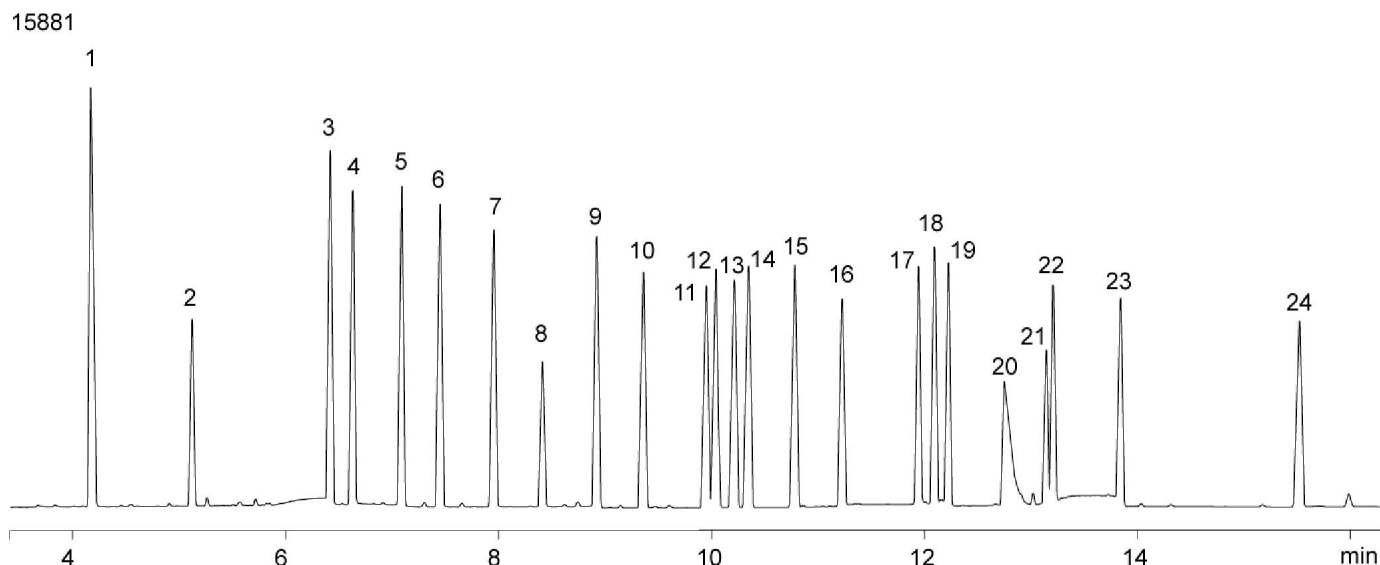Zebron™**Analyst Note:** All analytes are 20ppm

ANALYTES:

- | | | |
|---------------------------|-----------------------|-----------------------|
| 1 1-Bromo-2-Nitrobenzene | 10 Heptachlor Epoxide | 19 DDT |
| 2 Tetrachloro-m-xylene | 11 Endosulfan I | 20 Endrin Aldehyde |
| 3 a-BHC | 12 g-Chlordane | 21 Methoxychlor |
| 4 Pentachloronitrobenzene | 13 a-Chlordane | 22 Endosulfan Sulfate |
| 5 g-BHC (Lindane) | 14 DDE | 23 Endrin Ketone |
| 6 Heptachlor | 15 Dieldrin | 24 Decachlorobiphenyl |
| 7 Aldrin | 16 Endrin | |
| 8 b-BHC | 17 DDD | |
| 9 g-BHC | 18 Endosulfan II | |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.com Phenomenex products are available worldwide. Email us at international@phenomenex.com					
Australia tel.: 02-9428-6444 email: info@phenomenex.com.au	Austria 01-319-1301 anfrage@phenomenex.com	Canada (800) 543-3681 info@phenomenex.com	Denmark 4824 8048 dkinfo@phenomenex.com	France 01 30 09 21 10 franceinfo@phenomenex.com	Germany 06021-58830-0 anfrage@phenomenex.com
Ireland tel.: 01 247 5405 email: eireinfo@phenomenex.com	Italy 051 736176 italiainfo@phenomenex.com	New Zealand 09-4780951 info@phenomenex.co.nz	Puerto Rico (800) 541-HPLC info@phenomenex.com	United Kingdom 01625-501367 ukinfo@phenomenex.com	USA (310) 212-0555 info@phenomenex.com