

Tocopherols on Hydro-RP

Column: Synergi 4u Hydro-RP 80A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4375-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 80

Eluent B: Methanol


Pct. B: 20

Flow Rate: 3 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

Page 1 of 1

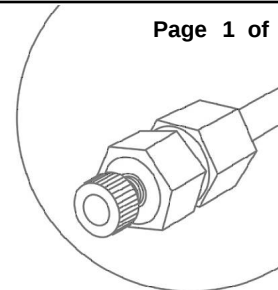

Products used in this Application:

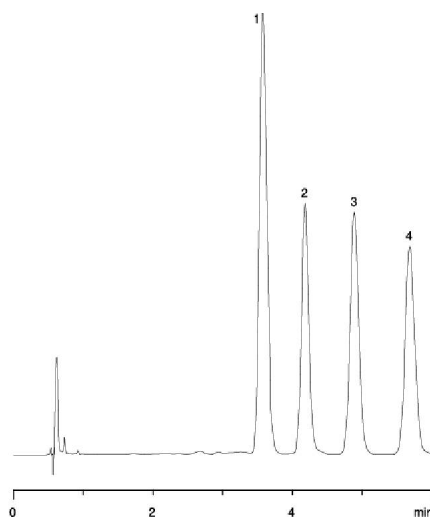synergi[™]

ANALYTES:

- 1 Delta tocopherol
- 2 Beta tocopherol + Gamma tocopherol
- 3 Alpha tocopherol
- 4 Tocopherol acetate



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.