


Methoxamine


Column: Chirex 3020 (S)-tert-leucine and (R)-1-(α -naphthyl)ethylamine
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3020-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (60:35:5) ethanol-TFA 20:1
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 290 nm (ambient)

Pct. A: 100

Page 1 of 1

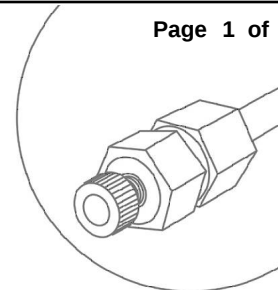

Products used in this Application:

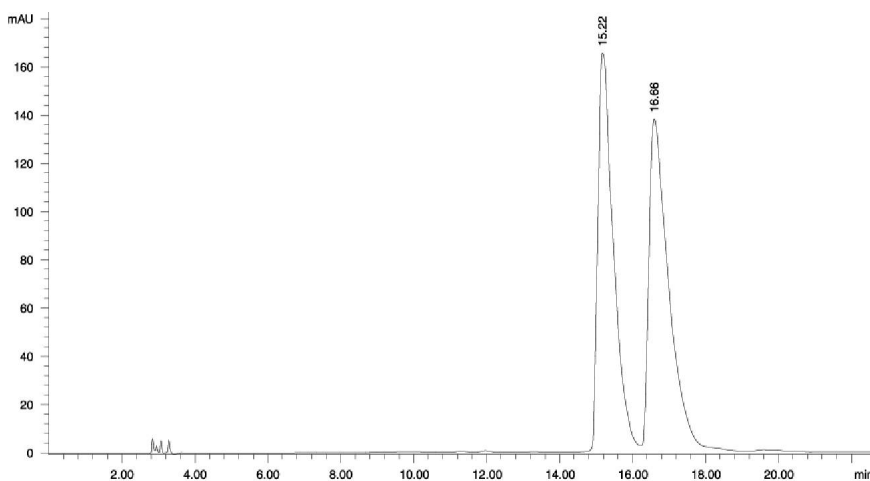CHIREX® $\alpha = 1.11$

ANALYTES:

- 1 Methoxamine
- 2 Methoxamine