

Complex Sugars on Luna Amino

Column: Luna 5u NH2 100A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4378-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Flow Rate: 3 mL/min


Col. Temp.: 40 °C

Detection: Refractive Index (RI) (40 °C)

Pct. A: 65


Pct. B: 35

Page 1 of 1

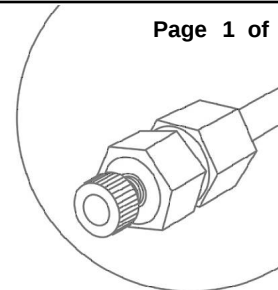

Products used in this Application:

LUNA®

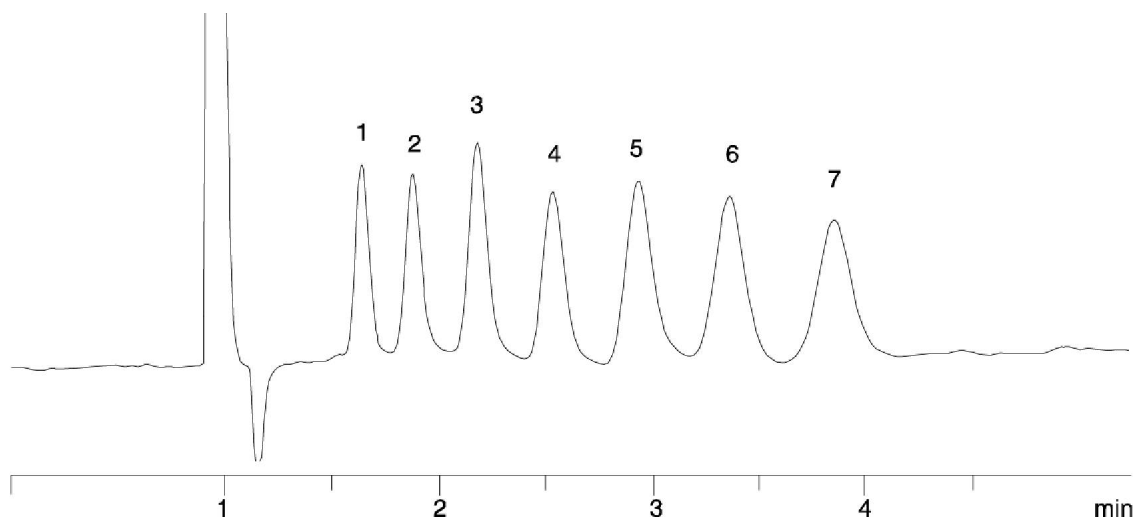

ANALYTES:

- 1 Glucose
- 2 Maltose
- 3 Maltotriose
- 4 Maltotetraose
- 5 Maltopentaose
- 6 Maltohexaose
- 7 Maltoheptaose