


Panax Ginseng on Synergi Max-RP

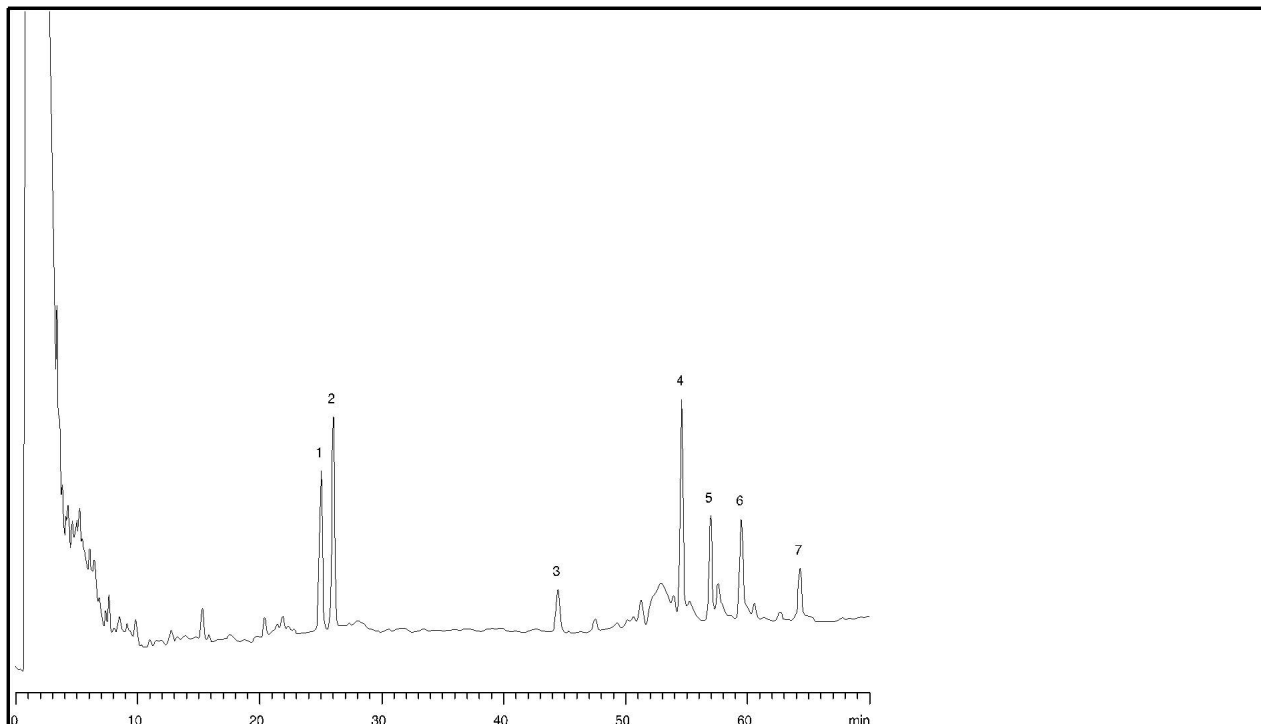

Column: Synergi 4 μ m C12 Max-RP

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	85	15
2	70	65	35

Flow Rate: 1.5 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 205 nm

SPE:

SAMPLE ANALYTE

- 1 *Ginsenoside Rg1*
- 2 *Ginsenoside Re*
- 3 *Ginsenoside Rf*
- 4 *Ginsenoside Rb1*
- 5 *Ginsenoside Rc*
- 6 *Ginsenoside Rb2*
- 7 *Ginsenoside Rd*

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14227