

EPA Method 632, Carbamate and Urea Pesticides

Page 1 of 1

Column: Luna 5u C18(2) 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

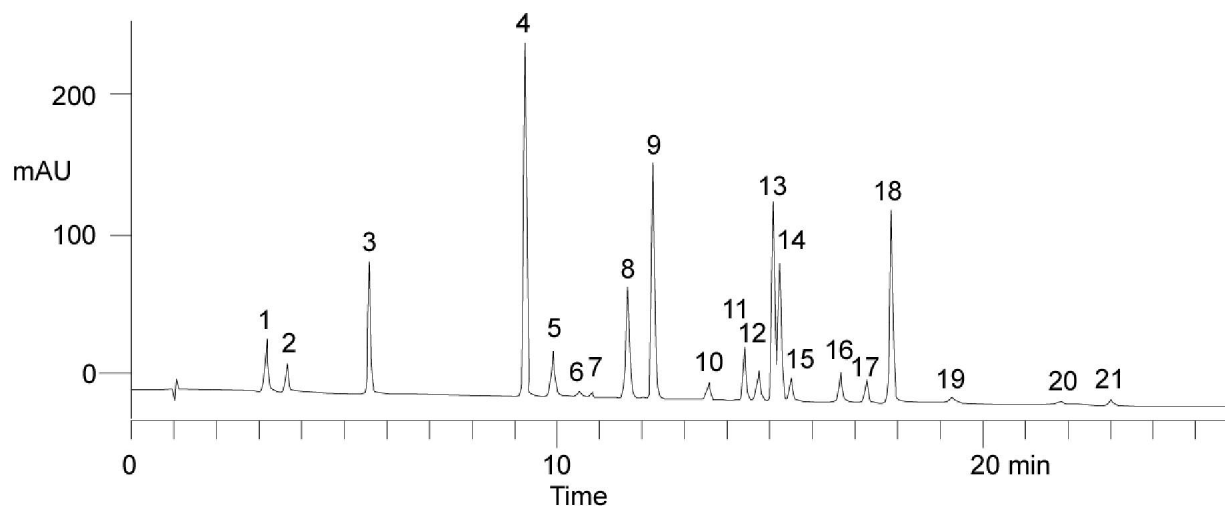
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	45	55
	2	6	45	55
	3	9	45	55

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Products used in this Application:

ANALYTES:

1 Aminocarb	14 Monuron
2 Barban	15 Monuron-TCA
3 Carbaryl	16 Neburon
4 Carbofuran	17 Oxamyl
5 Chlorpropham	18 Propham
6 Diuron	19 Propoxur
7 Fenuron	20 Siduron
8 Fenuron-TCA	21 Swep
9 Fluometuron	
10 Linuron	
11 Methiocarb	
12 Methomyl	
13 Mexacarbate	