

Carbamates on Envirosep-CM (EPA Method 632)

Page 1 of 1

Column: EnviroSep-CM

Dimensions: 175 x 3.2 mm ID

Order No: 00V-3133-R0

Elution Type: Gradient

Eluent A: Water



Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	14	50	50

Flow Rate: 1 mL/min

Col. Temp.: 45 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

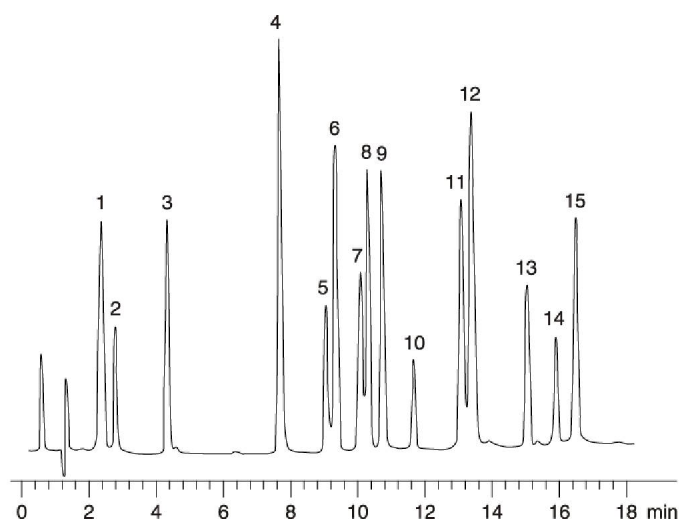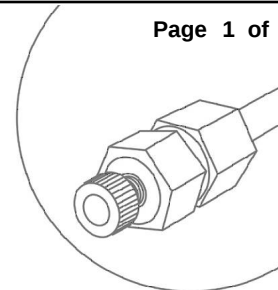

ANALYTES:

- | | |
|-----------------|------------|
| 1 Oxamyl | 15 Neburon |
| 2 Methomyl | |
| 3 Fenuron | |
| 4 Monuron | |
| 5 Propoxur | |
| 6 Carbofuran | |
| 7 Carbaryl | |
| 8 Fluometuron | |
| 9 Diuron | |
| 10 Protham | |
| 11 Siduron | |
| 12 Linuron | |
| 13 Chlorpropham | |
| 14 Barban | |