

Protein Standards

Column: Jupiter 5u C18 300A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4053-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in water


Eluent B: 0.1% TFA in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	20	35	65


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm (ambient)

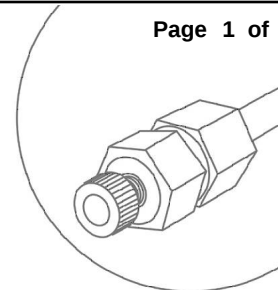

Products used in this Application:

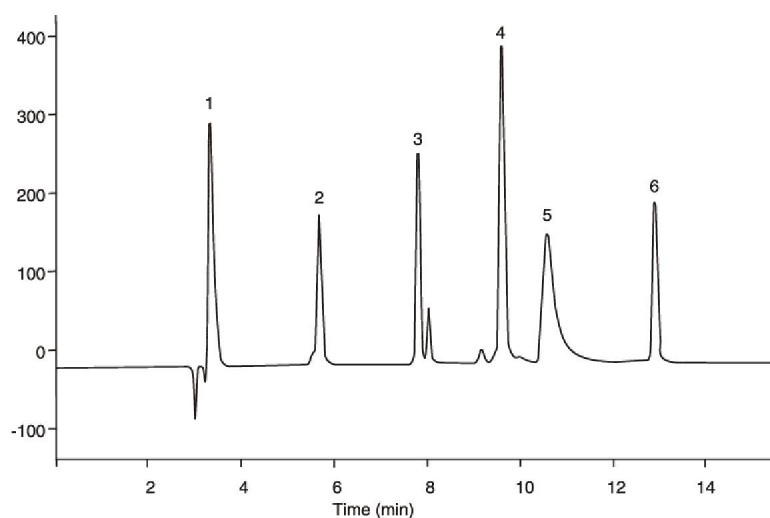

ANALYTES:

- 1 RNase
- 2 Insulin
- 3 An insulin contaminant
- 4 Lysozyme
- 5 Bovine serum albumin
- 6 Myoglobin