



Dextrans

Column: BioSep-SEC-S 4000**Dimensions:** 300 x 7.8 mm ID**Order No:** 00H-2147-K0**Elution Type:** Isocratic**Eluent A:** Water**Pct. A:** 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index (RI) (0 °C)

Page 1 of 1

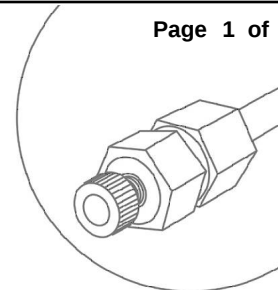

Products used in this Application:

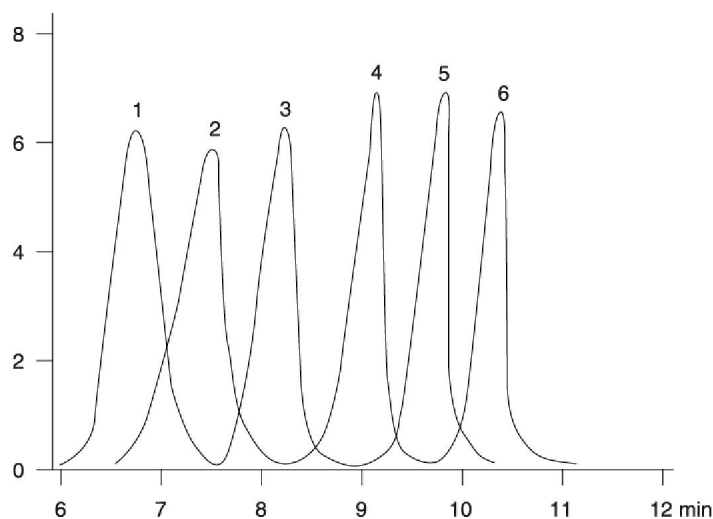Bio.Sep™**ANALYTES:**

- 1 Blue dextran
- 2 500,000
- 3 110,000
- 4 70,000
- 5 40,000
- 6 6,000