


농약 다성분 HPLC 분석

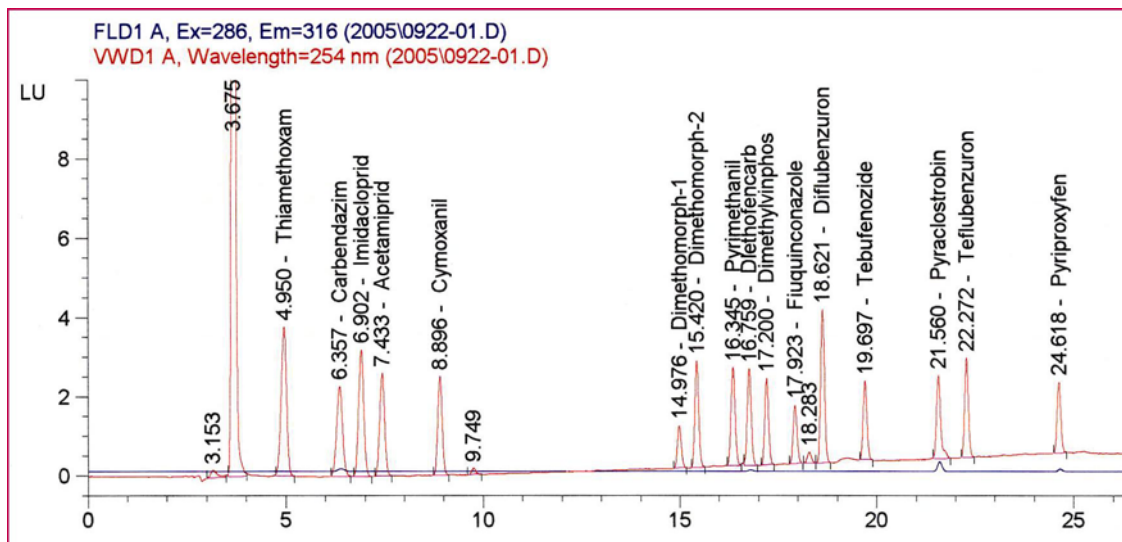

Chromatography Condition

Column : Gemini C18 5μ 250 x 4.6mm (Part No. 00G-4435-E0)
 Elution Type : Gradient
 Eluent A : Acetonitrile
 Eluent B : Water
 Detection : FLD Ex=286nm Em=316nm
 : VWD Wavelength=254nm
 Injection Volume : 10μL (Concentration 10 ppm)
 Temperature : 40℃

Gradient Profile	Time(min)	Flow rate(mL)	Solvent A	Solvent B
	0	1	30	70
	2	1	30	70
	22	1.2	85	15
	28	1.2	85	15
	30	1	30	70