


Peppermint Oil on ZB5ms 20x0.18x0.18

Page 1 of 1


Column: Zebron ZB-5MS
Phase: 5% Polysilarylene - 95% Polydimethylsiloxane copolymer
Dimensions: 20 meters x 0.18 mm x 0.18 µm
Order No: 7FD-G010-08
Oven Profile: 50°C to 100°C at 8°C/min to 225°C at 18°C/min for 4 min
Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: Split 100:1 1 µL @ 180 °C
Detection: Mass Selective (MSD) (180 °C)

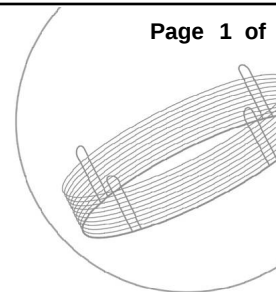

Products used in this Application:

Zebron™

15564

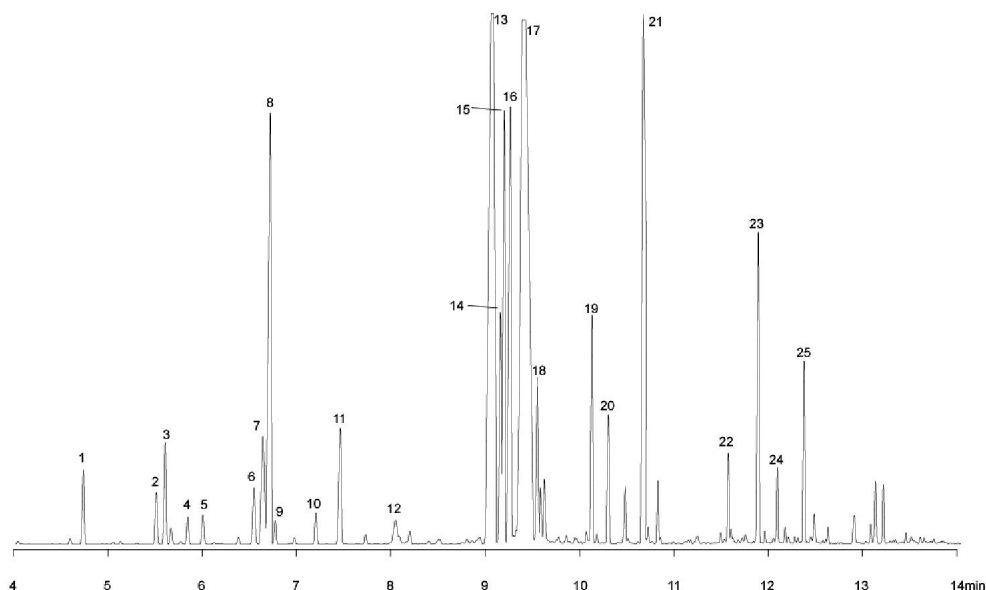

ANALYTES:

- | | | |
|------------------|-------------------|--------------------|
| 1 alpha-Pinene | 10 a-Terpinene | 19 Pulegone |
| 2 b-Phellandrine | 11 beta-Terpineol | 20 Piperitone |
| 3 beta-Pinene | 12 Linalool | 21 Menthyl Acetate |
| 4 b-Myrcene | 13 Isomenthone | 22 Bourbonene |
| 5 3-Octanol | 14 Menthonefurane | 23 Caryophyllene |
| 6 Cymene | 15 Menthone | 24 Farnasene |
| 7 D-Limonene | 16 Neo-menthol | 25 Germacrene |
| 8 Eucalyptol | 17 Menthol | |
| 9 trans-Ocimene | 18 Neoisomenthol | |

©Copyright 2000-2007 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555

www.phenomenex.com Phenomenex products are available worldwide. Email us at international@phenomenex.com					
Australia tel.: 02-9428-6444 email: info@phenomenex.com.au	Austria 01-319-1301 anfrage@phenomenex.com	Canada (800) 543-3681 info@phenomenex.com	Denmark 4824 8048 dkinfo@phenomenex.com	France 01 30 09 21 10 franceinfo@phenomenex.com	Germany 06021-58830-0 anfrage@phenomenex.com
Ireland tel.: 01 247 5405 email: eireinfo@phenomenex.com	Italy 051 736176 italiainfo@phenomenex.com	New Zealand 09-4780951 info@phenomenex.co.nz	Puerto Rico (800) 541-HPLC info@phenomenex.com	United Kingdom 01625-501367 ukinfo@phenomenex.com	USA (310) 212-0555 info@phenomenex.com